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Dear Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for the applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage and air.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Best regards,
Your Beurer Team

1. Getting to know your instrument

Function of the unit

This digital scale allows you to weigh yourself and to calculate your personal body mass index (BMI).

The scale is equipped with the following functions, which can be used by up to ten different people:

- Measurement of your body weight and display of time/temperature,
- Calculation of your average weight and your body mass index (BMI).

The memory function allows up to 10 people access up to 5 previously measured weight values.

2. Safety notes

Warning


- Do not stand on the outer edge to one side of the scale otherwise it may tip!
- Keep children away from packaging materials (risk of suffocation).
- The surface of the scale is slippery when damp or wet.

General notes

- The unit is for personal use only and is not intended for medical or commercial applications.
- Please note that technically related measuring tolerances are possible, as the scale is not calibrated for professional medical applications.
- If possible, always weigh yourself at the same time of day (preferably in the morning), after using the toilet, on an empty stomach and unclothed, in order to achieve comparable results.
- Load capacity: max. 180 kg (396 lb, 28 st), graduation for weight measurement: 100 g (0.2 lb).
- Repairs must only be carried out by Customer Services or authorised suppliers.
- Should you have any questions concerning the use of our devices, please contact your dealer or customer service.

3. Unit description

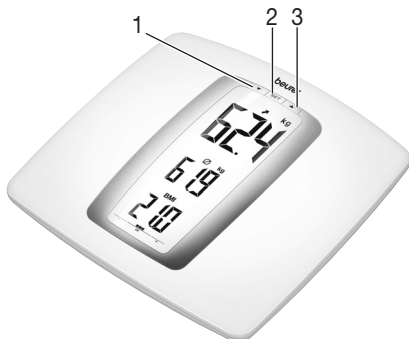- 1 Down button ▼
- 2 SET button
- 3 Up button ▲

4. Initial use

Changing the weight unit

In basic set-up, the scale displays values in kg. You can set the device to measure in kilograms (“kg”), pounds (“lb”) and stones (“st”) using the switch on the back of the device.

Inserting the batteries (3 x 1.5 V AA Micro)

Please remove the battery insulating strips and/or position the batteries correctly with regards to the polarity displayed in the battery compartment. If the scale fails to operate, remove the batteries completely and reinsert them.

Set time and date

The time can be changed by removing and reinserting the batteries or, if the scale is off, by holding the SET button for 3 seconds.

To set the time, follow these steps: hour – minute.

The setting to be changed will flash each time.

- To change a setting, press the buttons ▼ and ▲.
- To save the setting and move onto the next one, or to close the setting, press the SET button.

Note

- 24h display in kg mode
- 12h display in lb/st mode

Setting user data

The scale comes equipped with 10 locations for storing user data and space for more than 5 measurement values per user. You and, for example, the members of your family can save and later access the following personal settings in the storage locations:

User data	Parameters
Memory position	1 to 10
Sex	male  female 
Body size	100 to 220 cm (39” to 87”)
Age	5 to 100 years

Press SET and select the desired memory space using the ▼ or ▲ buttons. Press SET to confirm.

Similarly, when entering data proceed as follows: use ▼ or ▲ buttons to select and use SET to confirm.

Note: Without entering this data, it is not possible to calculate the BMI.

Setting up the scale

Set up the scale on a level and solid surface. A solid surface is a prerequisite for a correct measurement.

5. Operation

a. Weighing only

Step onto the scale (Fig. 1). Stand still. The time and temperature are displayed and the current weight is measured (Fig. 2). The measurement is complete when the weight unit stops flashing.

The scale switches itself off after approx. 15 seconds.

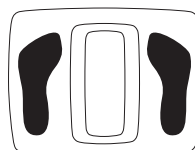
Fig. 1

Fig. 2

b. Weighing yourself and calculating BMI

- Using the ▼▲ buttons, select the user memory space where your settings have been saved.
- Wait until “0.0 kg” is displayed again. You can now use the scale. Stand still whilst weighing yourself.

The weight is displayed with trend, average weight and BMI with interpretation (Fig. 3).

The scale switches itself off after approx. 15 seconds.

Fig. 3

Accessing saved measurements

After a measurement with BMI calculation, your last 5 measurements can be displayed.

- After a BMI calculation, browse through the memory spaces for the measurements one by one using the ▼ button. Top: current weight; Middle: last weight; Bottom: memory space.

If no values are saved, no average value is displayed.

- By pressing ▼ again, the current weight, average weight and BMI are shown again.

6. Evaluation of results

Body mass index (BMI)

The body mass index (BMI) is a number that is often called upon to evaluate body weight. The number is calculated from body weight and height, the formula is – body mass index = body weight : height². The measurement unit for BMI is [kg/m²].

Man <20 years

- [Progress bar] +

Age	Under-weight	Normal weight	Over-weight
5	12,0–14,0	14,1–16,5	16,6–20,2
6	12,0–14,0	14,1–16,7	16,8–20,7
7	12,2–14,1	14,2–16,9	17,0–21,6
8	12,3–14,3	14,4–17,3	17,4–22,8
9	12,5–14,5	14,6–17,8	17,9–24,3
10	12,7–14,8	14,9–18,4	18,5–26,1
11	13,0–15,2	15,3–19,1	19,2–28,0
12	13,3–15,7	15,8–19,8	19,9–30,0
13	13,7–16,3	16,4–20,7	20,8–31,7
14	14,2–16,9	17,0–21,7	21,8–33,1
15	14,6–17,5	17,6–22,6	22,7–34,1
16	15,0–18,1	18,2–23,4	23,5–34,8
17	15,3–18,7	18,8–24,2	24,3–35,2
18	15,6–19,1	19,2–24,8	24,9–35,4
19	15,8–19,5	19,6–25,3	25,4–35,5

Woman <20 years

- [Progress bar] +

Age	Under-weight	Normal weight	Over-weight
5	12,0–13,8	13,9–16,8	16,9–21,3
6	12,0–13,8	13,9–16,9	17,0–22,1
7	12,2–13,8	13,9–17,2	17,3–23,3
8	12,3–14,0	14,1–17,6	17,7–24,8
9	12,5–14,3	14,4–18,2	18,3–26,5
10	12,7–14,7	14,8–18,9	19,0–28,4
11	13,0–15,2	15,3–19,8	19,9–30,2
12	13,3–15,9	16,0–20,7	20,8–31,9
13	13,7–16,5	16,6–21,7	21,8–33,4
14	14,2–17,1	17,2–22,6	22,7–34,7
15	14,6–17,7	17,8–23,4	23,5–35,5
16	15,0–18,1	18,2–24,0	24,1–36,1
17	15,3–18,3	18,4–24,4	24,5–36,3
18	15,6–18,5	18,6–24,7	24,8–36,3
19	15,8–18,6	18,7–24,9	25,0–36,2

Man/Woman ≥20 years

- [Progress bar] +

Age	Under-weight	Normal weight	Over-weight
≥20	15,9–18,4	18,5–24,9	25,0–40,0

Source: Bulletin of the World Health Organization 2007;85:660-7; Adapted from WHO 1995, WHO 2000 and WHO 2004.

- **Note that with very muscular physiques (body builders), the BMI interpretation will show overweight. The reason for this is that no account is taken of well above average muscle mass in the BMI formula.**
- **Remember that only long-term trends are important. Short-term fluctuations in weight over a few days are usually the result of a loss of fluid.**

Interpretation

The BMI values are interpreted in the evaluation graph as bars. When the values are low, the bar is shown to the left (-). When the bar is in the middle area, the values are normal. While in the right-hand area (+) the values are relatively high.

Trend display

This scale is equipped with a trend display for weight, which compares the current reading with the average value of all stored measurements.

- Upward trend display : The current reading is higher than the average of stored measurements.
- Downward trend display : The current reading is lower than the average of stored measurements.

7. Storage and maintenance

The precision of the measurements and service life of the device depends on its careful handling:

Important

- The unit should be cleaned occasionally. Use a slightly damp cloth to clean the device. You may also use a small quantity of cleaning agent if desired. Do not use abrasive detergents and never immerse the unit in water.
- Protect the scale from knocks, damp, dust, chemicals, drastic changes in temperature, electromagnetic fields and nearby sources of heat (ovens, heaters).

8. Disposal

Standard and rechargeable batteries should not be disposed of separately from the household waste. As a consumer, you are legally obliged to return used batteries for proper disposal. You can hand in your used batteries at public collection points in your district or sales outlets where batteries of this type are sold.

Note

The codes below are printed on batteries containing harmful substances:

Pb = Battery contains lead,

Cd = Battery contains cadmium,

Hg = Battery contains mercury.

For environmental reasons, do not dispose of the scale with the batteries in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point.

Observe the local regulations for material disposal.

Dispose of the unit in accordance with EC Directive **2002/96/EC – WEEE** (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

9. What if there are problems?

If the scale detects an error when measuring, the following is displayed.

Display	Cause	Remedy
OLd	The maximum load capacity of 180 kg was exceeded.	Do not load the unit with over 180 kg.
Lo	The batteries are almost empty.	Replace the batteries.

No measurement possible

Possible errors	Remedy
The batteries in the scale are empty.	Replace the batteries.

Incorrect weight measurement

Possible errors	Remedy
Scale is on carpeting.	Place the scale on a level and solid surface.
Incorrect zero point of scale.	Place a scale on a secure, flat surface. Switch on the scale and wait until the scale automatically switches off again (do not step onto the scale before then). Then switch on the scale again and repeat the measurement.

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