# beurer

## FM 35

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
**Dear customer,**

We are happy that you have decided on a product from our range. Our name stands for high-quality and exhaustively tested products from the areas of warmth, gentle therapy, blood pressure, body temperature, weight, massage and air.

Please read these instructions for use carefully and follow the directions.

Yours faithfully,  
The Beurer Team

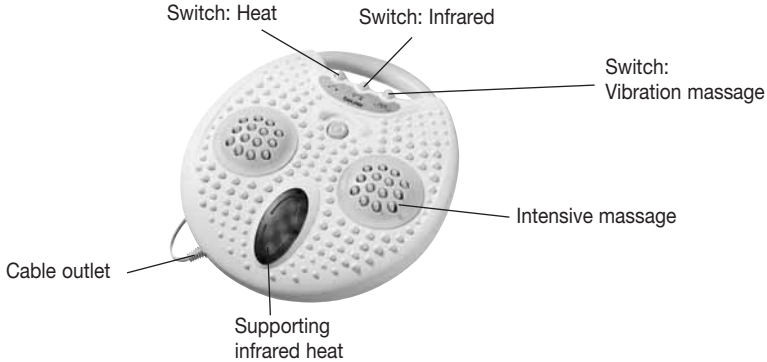**Application**

To avoid damage or premature aging of the foot massager, ensure that your feet are dry and clean when using it. Before operating the massager, ensure that the function switches for heat, infrared and massage are in central positions (position "0"). Place the massager in front of a chair, sit down and place your feet on the massager surface. Ensure that it is level. The heat, infrared and massager functions can be used together or separately

**Switch positions:**

<p><b>Heat:</b> 0 = switched off 1 = switched on</p>	<p><b>Infrared light:</b> 0 = switched off 1 = switched on</p>	<p><b>Massage:</b> 0 = switched off 1 = medium vibration massage 2 = maximum vibration massage</p>
--	--	--

We recommend a session of not more than 20 minutes. During the vibration massage, the heated surface and supporting infrared heat additionally relax your tense foot musculature and stimulate blood circulation. Because there may be a continuous tingling sensation at the soles of the feet after a massage, we advise you to begin with the medium vibration massage (switch position 1).

**If the surface is too warm for you, switch off the massager for several minutes to allow it to cool.**  
**Please note that massage is not a substitute for medical treatment.**  
**In the case of inflammation, chronic symptoms or pain, check with your GP before using the massager.**

**Cleaning**

Always remove the power plug from the socket outlet before cleaning. Use a damp cloth to clean the plastic case. You can also clean the surface with mild washing-up liquid. Never use abrasive, caustic or solvent-containing cleaners. When cleaning, ensure that no moisture enters the inside of the case.

**Important information – please keep for future reference**

- The device is intended to be used exclusively for the purpose described in the manual.

- The manufacturer cannot accept any liability for damage caused by improper use.
- Check the device for signs of damage or wear before use.
- In case of damage or faults, switch off the device immediately and unplug from the mains.
- If the mains cable of the device is damaged, it must be replaced by the manufacturer or an authorised aftersales service centre.
- Unprofessional repairs can pose significant risks to the user.
- Never leave the device on in your absence.
- Do not stand on the device.
- Never use the device in the bath, shower or wash basin and never immerse it in water.
- Children are unaware of the potential risks associated with electrical appliances. You should therefore make sure that the device cannot be used by children without supervision.
- The massager should not be operated continuously for longer than 20 minutes.
- Consult your doctor if you have any medical condition which affects the legs or feet, e.g. phlebitis, varicose veins or diabetes.
- The surface of the appliance is hot, persons who are insensitive to heat must use it with particular care.
- Do not open the foot massager. Failure to comply invalidates the warranty. You must not adjust or repair the foot massager yourself or its correct function will no longer be guaranteed.

Should you have any further queries regarding the application of the device, please contact our after-sales service staff, who will be glad to advise you.

## Disposal

Please dispose of the device in accordance with the directive 2002/96/EG – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please refer to the local authorities responsible for waste disposal.

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