

Caremaxx®

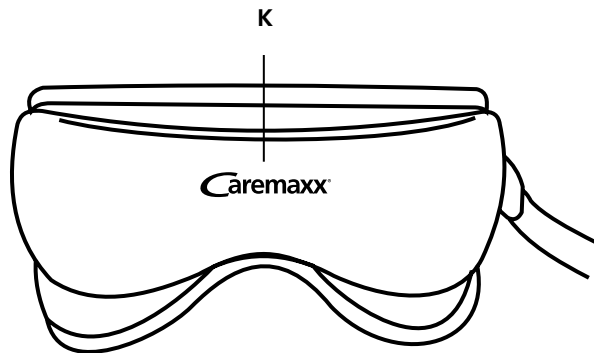- Ontspannende massage
- Massage de détente
- Relaxing massage
- Entspannende Massage

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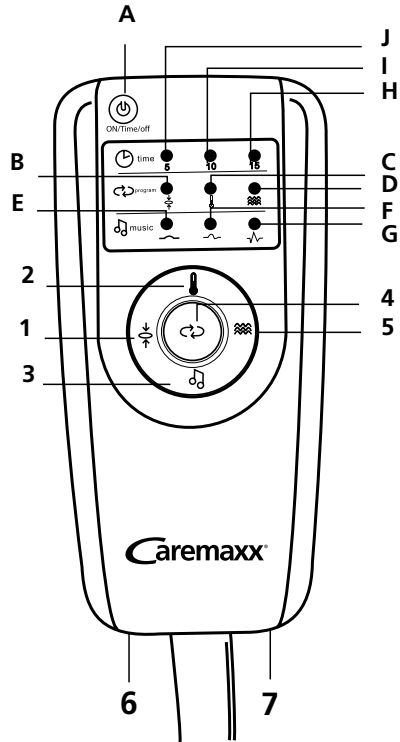E-mail : info@caremaxx.nl
Site Internet : www.caremaxx.nl

1 Caremaxx

1.1 Thank you for purchasing Caremaxx eye massager.

Carefully read these user instructions before using the eye massager.

1.2 Safety measures

- Only use this eye massager for purposes described in this user manual.
- The guarantee becomes void in the event of improper use.
- This device is not intended for commercial use.
- Do not leave this device unsupervised within the reach of children.
- Do not use the eye massager while talking, driving a car or while taking a bath or shower.
- Do not use the eye massager if you suffer from eye problems such as staring, retina infections etc.
- Repairs may only be carried out at authorised service points.
- Do not expose the device to extremely high or low temperatures or direct sunlight.
- Do not dip any part of the eye massager in water or any other liquid.
- Close your eyes during use. Keep them closed for a few moments after use, before opening them slowly.

Delivery and packaging

We want our quality products to reach you in the same condition that we sent them. That is why your **Caremaxx** eye massager comes in strong protective packaging which is made from environment-friendly materials. If you

do not want to use the box to store the device, then you can simply dispose of it like normal paper. If you discover transport damage while unpacking, then please contact your supplier.

The package contains following components:

- 1x eye massager with remote control
- 1x headphones
- 1x adaptor
- 4x AA batteries
- 1x user manual

2 Using your eye massager

2.1 Placing and removing the batteries.

Turn the remote control upside-down. Push the battery cover in the direction of the arrows in order to open it.

As shown, place 4 AA batteries with the plus (+) facing upwards and minus (-) facing downwards. Replace the cover until a click is heard.

In order to prevent leakage, remove the batteries if you do not use the eye massager for long periods.

2.2 Instructions for use

- Ensure that batteries have been placed correctly. It is also possible to use a mains adaptor. First remove the batteries and then plug the mains adaptor into the adaptor socket ⑥.
- Connect the headphones to the headphones input ⑦.
- Put on the eye massager and ensure that the elastic is the correct length, so that it is not too tight.

2.3 The massage

Turn on the eye massager by pressing the On/Off switch (A). The on/off indicator on the glasses will light up (K).

The massager starts with a 5-minute varied massage programme. This involves three variations: Air pressure + heat, vibration + heat, air pressure + vibration.

The device is standardly programmed with relaxing materials.

You can select the desired feature combination using the massage button ④.

Programme 1: Interchanging air pressure + heat, vibration + heat, air pressure + vibration.

Programme 2: Air pressure + heat

Programme 3: vibration + heat

Programme 4: Interchanging air pressure and vibration

Programme 5: Air pressure

Programme 6: heat

Programme 7: vibration

You can configure the air pressure massage in 3 ways by using the air pressure button ①:

1. Eye and temple massage
2. Eye massage
3. Temple massage

You can configure temperature in 3 ways by using the heat button ②:

1. Low temperature
2. High temperature
3. Off

You can configure the vibration in 10 ways by using the vibration button ⑤:

1. Strong vibration
2. Average vibration
3. Soft vibration
4. Interchanging average and soft vibration
5. Interchanging strong and soft vibration
6. Vibration interval - long
7. Vibration interval - medium
8. Vibration interval - short
9. Swing – interchanging left and right
10. Fast swing – interchanging left and right

2.4 Sound

Select your preferred sound using the music button ③.

You can choose from the following natural sounds:

1. relaxing music
2. trickling water
3. birds singing
4. the sea
5. crickets
6. falling droplets

Hold down the music button ③ for a few seconds in order to turn off the sound.

3 General

3.1 Cleaning and maintenance

Switch off the eye massager before cleaning them. Clean it using a dry cloth. Store the eye massager in a dry place, away from direct sunlight.

Remove the batteries from the device if it will not be used for long periods. This prevents the batteries from leaking.

3.2 Waste management

This instrument is not to be disposed of with the household waste. Consumers are obliged to send all electric or electronic instruments, whether or not they contain any hazardous materials, to a depository for environmental pollutants in their home town or to their specialized dealership to ensure their environmentally sound disposal. Do not forget to remove the battery from the instrument before throwing it away. Used batteries do not belong into the household waste, but should be taken to the pertinent waste management facility or a special collecting point at the supermarket or the electronics store. If you have any questions regarding waste management, please contact your municipality or your specialized vendor.

3.3 Warranty and repair conditions

If you wish to submit a warranty claim, please contact your supplier or our customer service. If the instrument is to be sent to our customer service, make sure to enclose a description of the complaint as well as a copy of your receipt of purchase. Warranty is granted subject to the following conditions:

1. The **Caremaxx** products are subject to a warranty period of three years, starting with the purchasing date as evidenced by the receipt or invoice.
2. All complaints due to material or manufacturing defects will be remedied free of charge within the warranty period.
3. A warranty case does not automatically imply an extension of the warranty period; neither for the instrument itself nor for the replaced parts.
4. Excluded from the warranty are:
 - a. Any damages caused by inappropriate use; for instance due to noncompliance with the user manual.
 - b. Damages caused by repair jobs performed by the purchaser or by any other unauthorized person.
 - c. Transport damages in transit accrued enroute from the seller to the consumer or from the consumer to the customer service.
 - d. Accessories that are subject to wear.
5. The producer accepts no liability for direct or indirect consequential damages caused by the instrument; not even if the damages to the instrument constitute an acknowledged warranty case.

3.4 Technical data

Product	Caremaxx eye massager
Model	30500
Vibration frequency	40Hz – 100Hz
Heating	30°C -45°C
Operating conditions	5°C - 40°C, <80% air moisture
Net voltage	AC: 220V +/-22V, 50Hz +/-1Hz DC: 6V
Remote control dimensions:	14.1 x 6.85 x 4.4 cm
Dimensions glasses:	23 x 10.5 x 10 cm
Battery	4x AA, 1.5V
EAN Code	87 17964 30500 4

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