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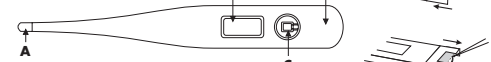
**HERSTELLER** Kangfu Medical Equipment Factory, 380 Ningqiang East Road, 325600 Yueqing, Zhejiang, China

# promed EN I Manual digital thermometer PFT-3.7

- A. Measuring tip
- B. Display
- C. Button
- D. Battery compartment cover
- E. Battery compartment
- F. Battery

**TECHNICAL DATA**  
**Device type:** Digital thermometer for oral, axillary and rectal measurements  
**Measuring range:** 32.0 °C - 42.0 °C  
**Scaling:** 0.1 °C  
**Measurement accuracy:** ± 0.1 °C at 35.5 °C - 42.9 °C  
± 0.2 °C below 35.5 °C and above 42.0 °C (at a room temperature of 25 °C)  
**Battery:** 1.5 VDC (LR41)

**Service life of the battery in continuous operation:** approx. 2,500 measurements  
**Automatic switch-off:** After 8 minutes  
**Class of protection:** IP 22  
**Dimensions in mm (LxWxH):** 128 x 19 x 10  
**Weight:** 10 g (incl. battery)

### DESCRIPTION AND APPLICATION OF THE DEVICE

Thank you very much for choosing the digital clinical thermometer **Promed PFT-3.7**. In order to ensure the quality and safety when using this device, this thermometer was designed and produced in accordance with Medical Devices Directive 93/42/EEC.

### SCOPE OF DELIVERY

- 1x Thermometer
- 1x Storage box
- 1x Manual

### BEFORE USING THE THERMOMETER FOR THE FIRST TIME

- Read these instructions of use with care and store them in a safe place.
- Make sure that the thermometer has not been damaged during transport.

### COMMISSIONING

- Prior to initial commissioning, please remove the contact protection strip of the battery. For this purpose, please proceed as described in section "Changing the Batteries".
- To switch the thermometer on, press the push button (C), all symbols will appear on the display (B) for approx. 2 seconds and you will hear a single beep.
  - The last measured value is now shown in the display for a moment (during this, a small M appears on the right edge of the display in addition). If the temperature is below 32 °C, you will then see a Lo and a flashing °C on the LCD display.
  - Position the sensor tip orally, axillary or rectally. Now, the thermometer will start the measurement automatically. The current temperature is continually shown during the measurement. As soon as the measurement is completed, you will hear 10 beeps and the measured temperature will be shown on the display. If the thermometer detects a temperature of > 37.7 °C, a high temperature alarm is given as 10 triple beeps.
  - Now you can switch off the thermometer by pressing the button (C) once. If you do not switch off the thermometer, it will automatically be switched off after 8 minutes.
  - The **Promed PFT-3.7** thermometer can be used for oral, axillary and rectal measurements.

### Oral measurement

- Lay the thermometer with the sensor tip into your tongue and close your mouth.
- Note: In order to avoid a inaccurate results, please keep your mouth closed during the measurement and breathe through your nose.
- The deviation from the actual body temperature may be -0.4 °C to -1.5 °C for oral measurement.

### Axillary measurement (measurement in the axilla)

- Clean your axilla before measurement. Lay the thermometer with the measuring tip in your axilla. Lower your arm and press it slightly towards your body during measurement.
- The deviation from the actual body temperature may be -0.5 °C to -2.0 °C for axillary measurement.
- Note: In order to prevent the measured values from being distorted, make sure to keep your armpit closed for 5 minutes before starting the measurement.

### Rectal measurement

- Although many patients perceive rectal measurement as unpleasant, it is the measuring method with no deviation to the actual body temperature.
- We recommend applying rect cream, lubricants or the like to the measuring tip before measurement.
- Lie on your side with slightly bent legs for measurement. If you use the thermometer for toddlers, we recommend laying the toddler on its stomach with its legs hanging down in an "L-position".
- Insert the thermometer approx. 2.3 mm into the rectum for measurement. Immediately stop if you sense any resistance when inserting the thermometer! Now you must measure until you hear the acoustic signal.
- If the thermometer has been used for rectal measurement once, it should not be used for any other measuring method due to hygienic reasons.

### MEMORY FUNCTION

The **Promed PFT-3.7** has a memory function that you can use to access the last measured value. Switch on the device by pressing the push button (C). All symbols of the display will appear for approx. two seconds and you will hear a short beep. Then, the last measured value is shown in the display for a moment (during this, a small M appears on the right edge of the display in addition).

### BATTERY CHANGE

If the battery power is too low, the symbol appears on the display (B). For battery change, you must remove the battery cover (D) from the battery compartment (E) and remove the old battery (F). Now you can insert a new battery (ensure correct polarity) into the battery compartment (E) and apply the battery compartment cover (D) to the thermometer.

### CLEANING AND DISINFECTION

The **Promed Thermometer PFT-3.7** is water-proof. Thus it is easy to clean and disinfect. Clean the area around the thermometer tip (approx. 5 cm from the measuring tip downwards) after each measurement with mild soap solution or with a cloth soaked in mild alcohol solution.

### SAFETY INSTRUCTIONS

- Never deposit or use the thermometer on hot surfaces, close to open gas flames or hot objects such as stove tops. Otherwise, the housing will melt.
- Use the thermometer only with batteries approved by the manufacturer.

- Use the thermometer indoors only.
- Keep the thermometer in a dry place.
- Do not bend the measuring tip.
- Do not bite the measuring tip during oral measurement.
- Some kind of supervision is required when using the thermometer around children.
- Make sure that children do not play with the thermometer. It includes small parts (battery and battery cover) which may cause suffocation when swallowed.
- For hygienic reasons, the thermometer should not be used for oral or axillary measurements after rectal measurement.
- Use the thermometer only for applications as specified in these instructions for use.
- The thermometer is exclusively intended for measuring human body temperature, as well as for the measuringspots in the human body as specified in the instructions for use.
- Using the thermometer in the vicinity of strong electromagnetic fields, e.g. next to a mobile telephone, may cause malfunctions.
- Please immediately contact your general practitioner in case of conspicuous results.
- This thermometer is not intended for use by persons (including children) with restricted physical, sensorial or mental abilities or lacking experience and/or lacking knowledge, unless they are supervised by a person responsible for their safety or if they were instructed on the use of the thermometer.
- Never allow children to play with the packing material, as there is a danger of suffocation.

### METROLOGICAL CHECK (FOR COMMERCIAL USE)

The thermometer was calibrated by the manufacturer for the duration of 5 years. For commercial use, a metrological check must be made at least every two years. The check is subject to a charge and can be accomplished by the responsible authorities or by authorised maintenance service agencies in compliance with the "Medical Device Operator Ordinance".

### DISPOSAL

Electric appliances, accessories and packaging should be recycled in an environmentally friendly manner.

#### Only for EU countries:

- Do not dispose of electric appliances in the household waste! According to the European Directive 2012/19/EU on Waste Electrical and Electronic Equipment and its transposition into national law, electric appliances which are no longer functional must be separately collected and recycled in an environmentally friendly manner. Within the EU, this symbol indicates that this product must not be disposed of in the household waste. Electronic waste contains recyclable materials which should be recycled in order to prevent damage to the environment and/or human health by uncontrolled waste disposal. Therefore, please dispose of old devices in suitable collection systems or return the device for disposal to the shop where you bought it. They will provide for adequate recycling.
- Rechargeable batteries:** Do not dispose of rechargeable batteries in the household waste or throw them into fire or water. Rechargeable batteries must be collected, recycled or disposed of in an environmentally friendly manner.
- Only for EU countries:** According to directive 91/157/EEC it is mandatory to recycle defective or used batteries. Batteries which can no longer be used may be returned to: Promed GmbH, Lindenweg 11, D-82490 Farchant.

### WARRANTY

This thermometer was manufactured with utmost care and checked before leaving the factory. Therefore this thermometer comes with a warranty of 24 months effective from the date of purchase and in accordance with the following terms and conditions.

- In case of verifiable material and/or manufacturing defects which may occur when using the thermometer in accordance with these instructions and which are identified within the warranty period, we will replace all defective parts of the thermometer free of charge within the warranty period including the share of labour costs for warranty repairs.
- Warranty will expire in case of the incorrect use of the thermometer (e.g. non-observance of the instructions of use), in case of interventions in the thermometer (e.g. opening the thermometer housing), as well as if spare parts are used which are not approved by **Promed**. Wear parts are not included in the warranty.
- The warranty period starts at the day of purchase. The filling of a warranty claim does not affect the warranty period. Warranty claims must be filed within the warranty period. Complaints after the expiration of the warranty period cannot be considered.
- Warranty claims within the scope of these warranty terms can only be filed if the date of purchase can be verified by a sales receipt or the like.
- Subject to technical and optical changes, as well as changes to the equipment!
- The statutory seller's warranty shall not be affected by our warranty terms.
- Promed** cannot be held liable for any possible translation errors.
- The following information is mandatory for trouble-free processing:
  - Original sales receipt or seller's stamp with date of purchase
  - Identified defect
  - Device designation / type

### EXPLANATION OF THE SYMBOLS

**Symbol for WARNING! This symbol refers to the safety-relevant chapters in these instructions of use!**

**Symbol for type B**

**Observe the symbol for instructions**

**Symbol for manufacturer and production site**

**Symbol for the date of manufacture**

CE 0123 **Symbol for the notified body**

**Symbol for the authorised representative in the European Union**

### TYPE PLATE

CE 0123

### CE DIRECTIVES

The requirements laid down in EU Directive 93/42/EEC in conjunction with the EU Directive 2007/47/EC for class IIa medical devices were met in accordance with standard IEC 60601-1:2012 +A1:2012. CE in accordance with the EU Directives 2014/35/EU, 2014/30/EU, 2011/65/EU, 1907/2006/EC, 2005/69/EC, 93/68/EEC and 2001/95/EC. The CE mark for the device refers to the EU Directive 93/42/EEC in conjunction with the EU Directive 2007/47/EC. Device classification: class IIa. Moreover, the device meets the requirements of the following standards: IEC 60601-1:2012 +A1:2012, EN 60601-2:2007+AC:2010, IEC 60601-1:11:2015, EN 124703

### MANUFACTURER

Kangfu Medical Equipment Factory, 380 Ningqiang East Road, 325600 Yueqing, Zhejiang, China

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