beurer

HK 62 Mobil

- Ⓛ **Akku-Heizkissen**
Gebrauchsanleitung
- Ⓜ **Battery-Operated Heating Pad**
Instruction for Use
- Ⓝ **Coussin chauffant à piles**
Mode d'emploi
- Ⓔ **Almohada eléctrica con acumulador**
Instrucciones para el uso
- Ⓜ **Termocuscino a batteria**
Istruzioni per l'uso
- Ⓝ **Akülü ısı yastığı**
Kullanma Talimatı
- Ⓜ **Акумуляторная электрогрелка**
Инструкция по применению
- Ⓜ **Poduszka elektryczna z akumulatorem**
Instrukcja obsługi

Inhaltsverzeichnis • Table of content • Sommaire • Índice • Sommario • İçindekiler • Оглавление • Spis treści

Ⓓ	Deutsch	Seite	4-9
ⒼⒷ	English	Page	10-15
Ⓕ	Français	Page	15-20
Ⓔ	Español	Página	20-26
Ⓘ	Italiano	Pagina	26-31
ⒹⓇ	Türkçe	Sayfa	32-37
ⓇⓊⓈ	РУССКИЙ	Стр	37-42
Ⓟ	Polski	Strona	42-47

Zeichenerklärung Etikett • Label legend • Légende de l'étiquette • Explicación de los símbolos de la etiqueta • Spiegazione dei simboli sull'etichetta • Şekillerin Açıklamaları Etiketi • Объяснение условных знаков на этикетке • Objasnienie oznaczeń etykiet

Hinweise lesen!
Read instructions!
Lire les instructions!
¡Leer las instrucciones!

Leggere le avvertenze!
Açıklamaları Okuyunuz!
Читайте указания по применению!
Czytać wskazówki!

Nicht gefaltet oder zusam-
geschoben gebrauchen!
Do not use when folded or creased!
Ne pas utiliser plié ou comprimé!
No utilizar estando plegado o
arrugado!
¡Non utilizzare ripiegato né accar-
tociato!

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Dear Customer

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage and air.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain. This diagnostic scale is designed to make a significant contribution to your health.

With kind regards
Your Beurer team

Items included in the package: 1 battery-operated heating pad, 1 battery pack, 1 charger, 1 car adapter, 1 belt extension, these instructions


1. Important safety information – read carefully and keep for later use

Failure to observe the information provided below can lead to personal injury or material damage (electric shock, burns, fire). The following safety and hazard information is not only intended to protect your health and the health of other persons but also to protect the product. Therefore, pay attention to this safety information and pass on these instructions together with the appliance when you sell it or pass it on.

- Only use this battery-operated heating pad for the purpose described in the instructions for use.
- This battery-operated heating pad is not intended for use in hospitals.
- This battery-operated heating pad is not intended for use by persons (including children) with restricted physical, sensory or mental capacity or with no experience and/or knowledge, unless they are supervised by a person who is responsible for their safety.
- It is not intended for use by helpless persons, small children or persons who are insensitive to heat (e.g. diabetics, persons with skin changes resulting from disease or scarring on areas of the skin in the area of application, after consuming pain-relieving medications or alcohol).
- Excessively long use can cause the skin to be burned.
- Caution! Never fall asleep while the battery-operated heating pad is switched on.
- Do not apply the battery-operated heating pad to parts of the body which are inflamed, injured or swollen. If in doubt, seek medical advice prior to use.
- The electric and magnetic fields emitted by this battery-operated heating pad may interfere with the operation of your pacemaker in some circumstances. However, they are far below the permissible limits: electric field strength: max. 5000 V/m, magnetic field strength: max. 80 A/m max., magnetic flux density: max. 0.1 millitesla. Please ask your doctor and the manufacturer of your pacemaker before using this battery-operated heating pad.
- The battery pack and connections must be protected against moisture, impacts or direct sunlight.

- The battery in the battery pack must not be disassembled, thrown into a fire or short-circuited.
- Only use the original battery charger (type P-050 B) to charge the battery pack. The battery charger is designed for indoor use only and must not be exposed to moisture.
- Use the battery-operated heating pad as follows
 - only connect it to the voltage specified on the appliance
 - only use it with the original AP2 battery pack supplied
 - do not operate it unattended
 - never switch on when folded or compressed
 - do not wedge it between objects
 - do not bend it sharply
 - do not use it with animals
 - do not use it when damp
- Do not operate it along with other heating devices such as heated car seats or heating underblankets.
- Never pull, twist or make any sharp kinks in the cables.
- Do not stick pins or pointed objects into this battery-operated heating pad.
- This battery-operated heating pad may not be used by children unless the switch has been preset by a parent or supervisory person, or the child has been sufficiently instructed on how to use this battery-operated heating pad safely.
- Children should be supervised to ensure that they do not use the heating pad as a toy.
- This battery-operated heating pad must be checked often for signs of wear or damage. If signs of this kind are discovered, or if this battery-operated heating pad has been used incorrectly, it must be checked by the manufacturer or dealer before it is used again.
- Repairs may only be carried out by specially trained personnel at a repair workshop authorized by the manufacturer, as special tools are required. Incorrectly conducted repairs can lead to considerable danger for the user.
- Should you still have questions on using our appliance, please contact our customer service.

2. Description of Device:

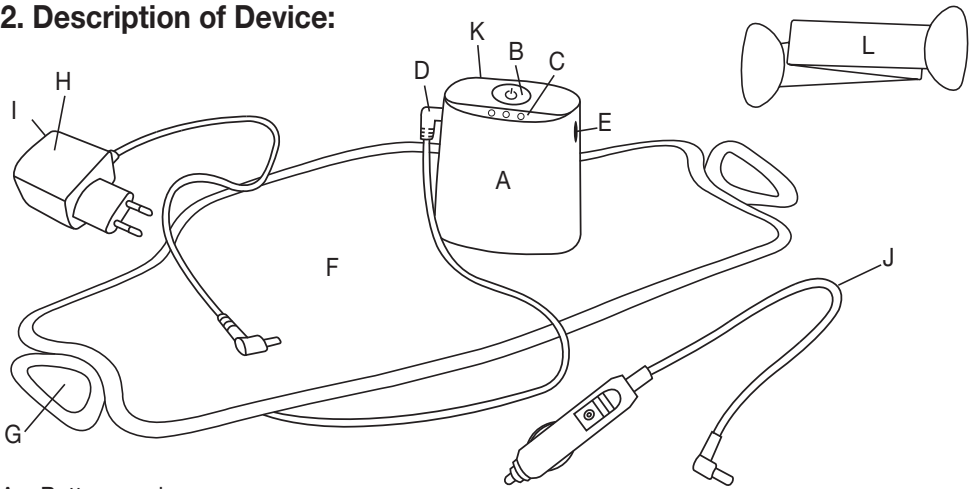- A Battery pack
- B Function key for temperature adjustment
- C Function display of battery pack
- D Connection for battery-operated heating pad
- E Connection for car adapter
- F Battery-operated heating pad
- G Elastic band with Velcro fastening
- H Charger
- I Function display of charger
- J Car adapter
- K Belt clip (removable)
- L Belt extension

3. Proper Use

This battery-operated heating pad is designed for warming the human body. You can use this heating pad to direct heat to specific parts of your body and you are no longer reliant on an electrical outlet.

The battery-operated heating pad provides sustained, soothing heat during any conceivable activity inside or outside.

This battery-operated heating pad is not intended for use in hospitals or for commercial applications. It must not be used for warming infants, small children, helpless persons or those insensitive to heat or for animals.

This battery-operated heating pad is made of breathable velour and is very adaptable and soothing to the skin.

You may perceive a plastic smell when using the battery-operated heating pad for the first time.

This will disappear after a short time.

This battery-operated heating pad can be used all purposes, such as on the stomach, back, joints and neck. If used for the neck area, the battery-operated heating pad can be rolled up.

The pull-out Velcro fastenings make it very easy to fasten it to the body. The elastic ensures that the battery-operated heating pad is always perfectly positioned.

4. Using it for the first time

4.1 Charging the battery pack

- Important: the battery pack must be charged prior to initial use or after prolonged periods without use.
- plug the charger into an electrical outlet with a voltage corresponding to the specifications on the charger. The green LED on the charger function display lights.
- Connect the charging plug to socket „E“ on the battery pack. The red LED „C“ on the function display of the battery pack lights red.
- Only charge the battery pack in a dry environment and at temperatures between 10 °C and 35 °C.
- The battery is fully charged when red LED „C“ on the function display of the battery pack is off (approx. 4 hours). This is followed by a maintenance charge of the battery pack.
- Please note that the full capacity of the battery pack is only reached after several charge and discharge cycles.
- The battery pack should be fully discharged and recharged at least twice a year to maintain its full capacity.
- Do not cover the battery pack with blankets, pillows or similar objects while it is charging. The battery pack will switch off if it is overheated.
- The heating pad cannot be used during charging.

4.2 Operating the battery-operated heating pad with the battery pack

- Remove the charger plug from the battery pack
- Connect the heating pad plug to connection D on the battery pack
- Select the required temperature setting with function key B.

4.3 Selecting temperature:

Press function key B to select the temperature.

The fastest way to heat up the heating pad is to set the switch to the highest setting first (level 2). You can switch back to a lower temperature setting later as required.

We recommend the lowest temperature setting for extended use of the heating pad.

4.4 Temperature settings:

- off = no LED on
- minimum heat = 1 blue LED on (press once)
- medium heat = 2 blue LEDs on (press twice)
- maximum heat = 3 blue LEDs on (press three times)

4.5 Automatic switch-off

The battery-operated heating pad has an automatic switch-off function. It stops the heat approx. 90 minutes after switching on the heating pad. Function key C switches off at this point.

4.6 Battery pack charge status

When the battery is discharged, red LED „C“ on the function display on the battery pack lights. The heating pad then no longer provides warmth. Charge the battery pack as described in 4.1. Please note that the battery pack should be protected against extreme cold (under clothing if applicable) to maintain its full efficiency.

If all three blue LEDs on the battery pack function display are flashing, the heating pad is not correctly connected to the battery pack or there is a fault. Check the connection. If the three blue LEDs continue to flash, please send the battery-operated heating pad and the battery pack for repair.

4.7 Using the battery-operated heating pad in a car (plugged into the 12-V cigarette lighter socket)


- Connect the heating pad plug to connection D on the battery pack
- Connect the car adapter to connection E on the battery pack.
- Connect the car adapter to the cigarette lighter socket in the car.
- Select a temperature setting (see 4.3/4.4)
- Please note that the battery-operated heating pad must not be used in conjunction with heated car seats.
- If the battery-operated heating pad is not used, it should be switched off (LED on the function display off) and the car adapter disconnected.
- The battery is not charged during car operation.

4.8 Replacing the battery in the battery pack

If the battery-operated heating pad only heats for a short time, we recommend replacing the

the battery with a new original spare part (see 4.9).

- Unscrew the back of the housing (4 Philips-head screws).
- Remove the battery from the housing.
- Disconnect the connection between the battery and the board by pulling the white plug connector apart.
- Connect the new battery to the board by connecting the two parts of the white plug connector.
- Insert the new battery into the housing.
- Make sure that the connector wire feeds through the opening in the board (otherwise it will be impossible to close the back of the housing completely).
- Close the back of the housing and screw in the 4 Philips-head screws tightly.
If you have any problems with replacing the battery, please contact our customer service.

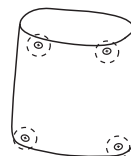4.9 Spare parts

You can obtain the following spare parts directly from Customer Service:

- Battery (Item No. 162.271 without housing)
- Battery charger (Item No. 546.006)
- Car adapter (Item No. 108.821)

5. Cleaning and care

Always disconnect the plug of the battery-operated heating pad from connector 0 of the battery pack before cleaning it.

Smaller stains on the battery-operated heating pad can be removed with a cloth or damp sponge and if necessary, some liquid detergent for delicate fabrics. Do not use cleaners containing solvents.

Please note that the battery-operated heating pad must not be dry-cleaned, wrung out, machine-dried, put through a mangle or ironed.

Set the washing machine at a gentle washing cycle at 40°C.

Machine-wash the battery-operated heating pad only in a cloth bag or pillow case.

For ecological reasons, we recommend that the heating pad is only washed together with other textiles. Use a detergent for delicate fabrics in quantities specified by the manufacturer.

Please note that excessive washing will wear out the heating pad. The heating pad should only be machine-washed a maximum of 10 times throughout its entire life.

Never use any clothes pegs or similar to secure the battery-operated heating pad when drying and never switch it on to dry it.

Do not connect the battery-operated heating pad to the battery pack until the plug and heating pad are completely dry.

6. Storage

If you do not use the battery-operated heating pad for longer periods, we recommend storing it in the original packaging in a dry place without weighting it down. Allow the battery-operated heating pad to cool down beforehand.

7. Disposal

Please dispose of the device in accordance with the Ordinance on Waste Electrical and Electronic Equipment (WEEE - 2002/96/EC). If you have any questions, please contact the municipal authorities responsible for waste disposal in your area. The battery must be removed from the battery pack first (see 4.5) and also properly disposed of.

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