# beurer

## Wärme-Unterbett TS 20, 23

TS 20

TS 23

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**Dear Customer,**

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for the applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage and air. Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Best regards,  
Your Beurer Team

**Items included in the package:** 1 Electric blanket, 1 switch, 1 operating instructions

# 1. Important information – keep for later use

Failure to observe the information provided below can lead to personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information not only serves to protect your health and the health of others, it also protects the product from damage. Therefore, observe this safety information and pass on these instructions together with the product.

- Only use this heated underblanket for the purpose described in these instructions.
- This heated underblanket is not intended for use in hospitals.
- This underblanket is not intended for use by persons (including children) with restricted physical, sensory or mental capabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety.
- Do not use this heated underblanket for helpless persons, small children or persons insensitive to heat.
- The electric and magnetic fields emitted by this electrical heated underblanket may possibly interfere with the operation of pacemakers. However, they are well below the permissible limits: electric field strength: maximum of 5,000 V/m, magnetic field strength: maximum of 80 A/m, magnetic flux density: maximum of 0.1 millitesla. Please therefore ask your doctor and the manufacturer of your pacemaker before using this product.
- This heated underblanket should
  - only be connected to the voltage specified on the appliance,
  - this heated underblanket should only be operated with the switch indicated on the device (label),
  - not be switched on if it is folded or compressed,
  - not be tucked into the bed or laid around the mattress,
  - never be wedged (especially with adjustable beds or divided mattresses),
  - never have any sharp creases,
  - never be used with animals,
  - never be used when damp.
- While this heated underblanket is switched on
  - never place any objects on it (e.g. suitcases or laundry basket),
  - never place any heat source on it, such as a hot water bottle, heating pad etc.

- When in use, the electronic components in the switch of the heated underblanket make the switch slightly warm. The switch should therefore never be covered or placed on the heated underblanket when it is in operation.
- Do not pull, twist or make any sharp kinks in the cables.
- Do not stick pins, needles or any sharp objects into this heated underblanket.
- This heated underblanket may not be used by children unless the switch has been preset by a parent or a supervisory person, or the child has been sufficiently instructed in how to use this heated underblanket safely.
- Children should be supervised to ensure that they do not use the device as a toy.
- This heated underblanket must be checked frequently for signs of wear or damage. If there are such signs or if the heated underblanket has been used incorrectly, it must be taken to either the manufacturer or the dealer before further use.
- Repairs may only be carried out by specially trained personnel or at an authorised repair workshop, as special tools are required. Incorrect repairs may cause considerable danger to the user.
- Should you still have questions on using our appliances, please contact our customer service.

## 2. Use

This electric underblanket is intended for warming beds. The electric underblanket heats up best when it is switched on when covered with a quilt. The heating up time with the electric underblanket covered is approx. 30 minutes. You can use either side of the electric underblanket.

The electric underblanket is fitted with the BEURER SAFETY SYSTEM (BSS). This electronic sensor technology prevents the blanket from overheating over its entire surface by switching it off automatically in the event of a fault. If, as a result of a fault, the electric underblanket is switched off automatically by the Beurer Safety System (BSS), the signal lamp will no longer light up even when the blanket is switched on. Please note that for safety reasons the electric underblanket must no longer be operated after a fault and must be sent in to the service address stated.

### 2.1 Operation

Place the electric underblanket on your mattress and spread it out flat starting from the foot end. Then as usual place the quilt over it, so that the electric underblanket is between mattress and quilt. Always ensure that the electric underblanket is spread out flat and is not folded or creased when in use. When making the bed, check its position – the electric underblanket must be switched on only when it is spread out flat. We urgently recommend that you switch on the electric underblanket for approx. 30 minutes before going to bed and covering it with a quilt to prevent the heat from escaping. The electric underblanket with removable switch must be operated only with the types of switch specified on the blanket.

**Attention: The electronic components in the switch will cause the switch to warm up slightly when the electric underblanket is being used. For this reason, the switch must never lie on the electric underblanket or be covered by blankets, pillows etc.**

### 2.2 Temperature selection

To operate the underblanket first connect the switch to the radiator by plugging in the connector. Plug the power cord into the mains outlet.

The quickest way to heat up the heated underblanket is to first select the highest setting. If it is to be used for several hours, we then recommend setting the lowest temperature.

## 2.3 Temperature levels TS 20, 23

0 = off

1 = minimum heat

2 = medium heat

3 = maximum heat


## 3. Cleaning and care

Always remove the power plug from the power supply before cleaning the electric underblanket. Small spots can be removed with a cloth or damp sponge and if necessary a little mild liquid detergent. Note that the electric underblanket must not be dry cleaned, wrung out, machine dried or ironed. **On no account switch on the electric underblanket to dry it.** Do not hang up the electric underblanket with clothes pegs or similar to dry it.

### 3.1 Cleaning and care of electric underblankets with removable switches

These electric underblankets are machine washable. Before the washing cycle, first disconnect the plug connection and thus remove the plug from the electric underblanket (see illustration).

Set the washing machine to a delicate washing cycle at 30 °C. Use a mild detergent with dosage according to the manufacturer's instructions. Only reconnect the switch to the electric underblanket after the plug connection and the electric underblanket are completely dry. Note that the electric underblanket will suffer from too frequent washing. For its entire life, therefore, the electric underblanket should be washed not more than 5 times in a washing machine.
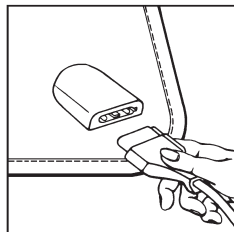
### 3.2 Protective covers

Beurer offers a white cotton protective cover (size 80 x 150 cm, Art. No. 160.100). A moisture-proof protective cover for incontinence is available in the size 80 x 150 cm (Art. No. 160.101). You can obtain protective covers from specialist dealers or directly from Beurer customer service.

### 3.3 Storage

If you are not using the electric underblanket for a long period of time, we recommend that you keep it in its original packaging in a dry place with no weight placed on it. Allow the underblanket to cool down beforehand.

### 3.4 Disposal

Please dispose of the blanket in accordance with the directive 2002/96/EG – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please refer to the local authorities responsible for waste disposal.

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