beurer FM 100

- (D) Wadenmassage-Gerät**
Gebrauchsanleitung
- (GB) Calf massager**
Instruction for Use
- (F) Appareil de massage des mollets**
Mode d'emploi
- (E) Aparato de masaje de pantorrillas**
Instrucciones para el uso
- (I) Massaggiatore per polpacci**
Istruzioni per l'uso
- (TR) Baldır masajı cihazı**
Kullanma Talimatı
- (RUS) Прибор для массажа икроножных мышц**
Инструкция по применению
- (PL) Masażer łydek**
Instrukcja obsługi

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Dear Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for the applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage and air.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Best regards,
Your Beurer Team

Items included in the package:

- Calf massager
- These operating instructions

Explanation of symbols

The following symbols are used in the operating instructions and on the type plate:

WARNING Warning of risks of injury or health hazards.

CAUTION Safety information about possible damage to appliance/accessories.

Note Important information.

The device has double protective insulation and therefore complies with protection class 2.

1. Getting to know your appliance

This device offers a powerful massage for the feet and calves.

The device is only for **private use** and is not intended for medical or commercial purposes.

The device offers the following advantages:

- 3 options for use: calf massage, foot massage and footrest
- 2 different massage functions in one: vibration and kneading massage

2. Proper use

This device is solely intended for massaging the calves and feet of human beings. It cannot replace medical treatment. Do not use the massage device if any of the following warnings apply to you.

WARNING

Use of the massage device is prohibited in the following cases:

- not with changes caused by disease or injury in the leg region (e.g. brittle bones, open wounds, varicose veins),
- not on children,
- not during pregnancy,
- never while you sleep,
- not on animals,
- not during activities during which an unexpected reaction may be dangerous,
- not while ingesting substances that reduce reaction times (e.g. painkiller medications, alcohol),
- never longer than 20 minutes (risk of overheating) and allow it to cool off for at least 20 minutes before using it again.

Please consult your doctor if you are unsure whether or not the massage device is suitable for you.

This device is not intended for use by persons (including children) with restricted physical, sensory or mental capabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety or they are given instructions on how to use the device.

Children should be supervised to ensure that they do not use the device as a toy.

You should consult your doctor before using the massage device if:

- if you suffer from a serious illness or are recovering from a leg operation,
- you are fitted with a cardiac pacemaker, implant or other aid,
- Thromboses,
- Diabetes,
- in case of pains of unclarified origin.

The device is only intended for the purpose described in these operating instructions. The manufacturer cannot be held liable for damage occurring as a result of inappropriate or careless use.

3. Safety information

Please read these operating instructions carefully. Failure to observe the information provided below can lead to personal injury or material damage. Keep these instructions in a safe place where they can also be read by other persons using the device. Always pass on these instructions together with the appliance if you sell it or pass it on.

WARNING

Keep children away from the packaging materials. Risk of suffocation.

WARNING: Electric shock

Just like any other electrical appliance, this massage device must be handled with due care and attention in order to prevent the risk of electric shock.

For this reason, the device must only be operated as follows:

- only with the supplied power supply unit at the mains voltage indicated on the power supply unit,
- do not use the device if damage to the device itself or any of its accessories is evident,
- never use the appliance during a thunderstorm.

In the event of a defect or malfunction, switch the device off immediately and disconnect it from the power supply. Never pull on the power cable or on the appliance in order to disconnect the plug from the mains outlet. Never hold or carry the device by the power cable. Keep the cables away from hot surfaces. Do not pinch, kink or twist the cables. Do not stick in any pins, needles or sharp objects. Never open the device for any reason.

Make sure that the massage device, the switch, the power supply unit and the cable do not come into contact with water, steam or other fluids.

Use the device only in dry inside rooms (e.g. never in the bathtub or sauna).

Never try to retrieve the device if it falls into water. Disconnect the plug immediately from the mains outlet.

Do not use the device if damage to the device itself or any of its accessories is evident.

Do not subject the device to any impact and do not drop it.

Repairs

WARNING

- Repairs to electrical appliances must only be carried out by qualified persons. Incorrect repairs can lead to considerable danger for the user. In the event of repairs, please contact our customer service or an authorized dealer.
- If the power cord is damaged it must be replaced by an authorised dealer.

Fire hazard

WARNING

There is a risk of fire if the device is used incorrectly or if it is not used according to these operating instructions.

For this reason, the massage device must only be operated as follows:

- Never use the device under covers, such as blankets or pillows.
- Never use the device in the vicinity of petrol or other easily flammable materials.

Handling

CAUTION

Switch the device off and disconnect the power supply after every use and before every cleaning.

- Make sure that all moving parts can move freely at all times.
- Never sit on the device or place your full weight on it. Do not place any objects on the device.
- Keep the device away from high temperatures.

Disposal



CAUTION

Follow local regulations regarding the disposal of materials. The appliance should be disposed of according to Regulation 2002/96/EC-WEEE (Waste Electrical and Electronic Equipment). In case of queries, contact the municipal authorities responsible for waste disposal in your area.


4. Appliance description

Calf massager

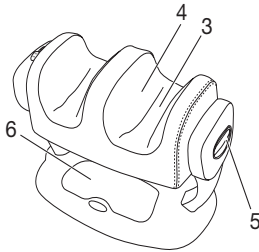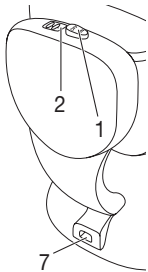
- 1 Button
- 2 Indicator light
- 3 Vibration massage
- 4 Kneading massage
- 5 Handle for setting the massage device to 3 different positions (see below)
- 6 Cord storage box
- 7 Socket for power cord

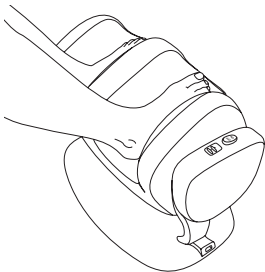Foot massage

Calf massage

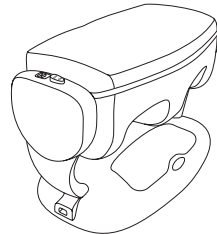Footrest

5. Operation

Connect the device to the mains by first plugging the power cord into the socket on the device and then plug it into the mains outlet.

The massage must be pleasant and relaxing at all times. Stop the massage or change your position or the contact pressure if you find it is painful or unpleasant.

- Place the device in front of the chair or couch and select one of the handle positions: footrest, foot massage or calf massage with the handle on the side. Sit comfortably on the seat cover. Switch on the device with the  button. Place your feet on the device.
- The device has an automatic shutoff. It will switch off automatically after 20 minutes of operation. Allow it to cool down for about 20 minutes before switching it on again. You should not exceed 30 minutes a day of massage. Using the device for longer can over-stimulate the muscles and cause tension instead of relaxation.
- After use always unplug the cord from the wall outlet, then disconnect it from the device and store it in the cord storage box on the device.
- DO NOT disconnect the device from the mains supply during operation. Always switch it off only once it has reached the starting position at the bottom

6. Care and storage

CLEANING

WARNING

- Disconnect the device before cleaning it.
- Only clean the appliance in the manner specified. Fluid must never be allowed to enter into the device or accessories.
- Smaller stains can be removed using a cloth or moist sponge together with a little liquid detergent for delicate fabrics if necessary. Do not use any cleaning agents containing solvents.
- Do not use the device again until it is completely dry.

Storage

If you do not use the device for longer periods, we recommend storing it in the original packaging in a dry place without weighting it down.

7. In the event of problems

Problem	Cause	Remedy
Massage unit moves more slowly.	The load on the massage unit is too great.	Reduce the pressure on your feet or calves.
Massage unit does not move.	Appliance not connected to the power supply.	Connect the plug to the mains supply and switch on the appliance.
	Device not switched on.	Switch the device on with  .
	Overheating protection tripped, device is damaged.	Contact your authorized dealer or the service address.

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