beurer MG 310

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Items included in the package:

- Shiatsu Massage Chair
- Removable covers and neck cushion
- Power supply unit
- These operating instructions

Explanation of symbols

The following symbols are used in the operating instructions and on the type plate:

WARNING Warning of risks of injury or health hazards.

CAUTION Safety information about possible damage to appliance/accessories.

Note Important information.

The device has double protective insulation and therefore complies with protection class 2.

Only use indoors.

1. About Shiatsu massage

Shiatsu is a form of body massage which was originally developed in Japan around 100 years ago. It is based on ideas anchored in traditional Chinese medicine, in particular the meridian system, i.e. the paths in the human body along which energy is conveyed.

The aim of any Shiatsu massage is to promote the physical, emotional and mental wellbeing of the person being treated. To do this, energy blocks and bottlenecks in the meridians need to be dissolved and the self-regulating ability of the body needs to be stimulated.

Shiatsu therapists achieve this by exercising pressure in flowing movements along the meridians. Although Shiatsu can be literally translated as “finger pressure”, practitioners also use the heels of their hands, their elbows and knees as well as their fingers. In the process, the therapist will use his or her full body weight in order to generate the required pressure.

2. Getting to know your appliance

Your Shiatsu massage device has rotating massage heads which imitate the pressure and kneading motions of a Shiatsu massage.

This device, with its vertically moving and rotating massage heads, offers a powerful, intensive deep massage over the entire back.

The device is only for **private use** and is not intended for medical or commercial purposes.

In addition, the Shiatsu massage device offers the following advantages:

- powerful Shiatsu massage along the spine with localised massage,
- rolling massage along the spine,
- adjustable massage zones (upper, lower or the entire back),
- adjustable massage intensity,
- automatic shutoff function after 15 minutes,
- with illuminated massage heads and adjustable heating function,
- removable and washable covers (seat, back, arm rests) and neck cushion,
- separate power supply unit.

3. Proper use

This device is solely intended for massaging the backs of human beings. It cannot replace medical treatment. Do not use the massage device if any of the following warnings apply to you.

WARNING

Use of the massage device is prohibited in the following cases:

- not with changes caused by disease or injury in the back region (e.g. disc problems, open wounds),
- not on children,
- not during pregnancy,
- never while you sleep,
- not in the car,
- not on animals,
- not during activities during which an unexpected reaction may be dangerous,
- not while ingesting substances that reduce reaction times (e.g. painkiller medications, alcohol),
- never longer than 15 minutes (risk of overheating) and allow it to cool off for at least 15 minutes before using it again.

This product can only be used with a maximum weight of 110 kg (250 lbs).

Please consult your doctor if you are unsure whether or not the massage device is suitable for you.

This device is not intended for use by persons (including children) with restricted physical, sensory or mental capabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety or they are given instructions on how to use the device.

Children should be supervised to ensure that they do not use the device as a toy.

The device has an adjustable heating function in the massage heads. Persons who do not feel heat must be careful when using this function.

You should consult your doctor before using the massage device if:

- you suffer from a serious illness or are recovering from an operation on the upper body,
- you are fitted with a cardiac pacemaker, implant or other aid,
- Thromboses,
- Diabetes,
- in case of pains of unclarified origin.

The device is only intended for the purpose described in these operating instructions. The manufacturer cannot be held liable for damage occurring as a result of inappropriate or careless use.

4. Safety information

Please read these operating instructions carefully. Failure to observe the information provided below can lead to personal injury or material damage. Keep these instructions in a safe place where they can also be read by other persons using the device. Always pass on these instructions together with the appliance if you sell it or pass it on.

WARNING

Keep children away from the packaging materials. Risk of suffocation.

WARNING: Electric shock

Just like any other electrical appliance, this massage device must be handled with due care and attention in order to prevent the risk of electric shock.

For this reason, the device must only be operated as follows:

- only with the supplied power supply unit at the mains voltage indicated on the power supply unit.
- Do not use the device if damage to the device itself or any of its accessories is evident.
- Never use the appliance during a thunderstorm.

In the event of a defect or malfunction, switch the device off immediately and disconnect it from the power supply. Never pull on the power cable or on the appliance in order to disconnect the power supply from the mains outlet. Never hold or carry the device by the power cable. Keep the cables away from hot surfaces. Do not pinch, kink or twist the cables. Do not stick in any pins, needles or sharp objects.

Make sure that the massage device, the switch, the power supply unit and the cable do not come into contact with water, damp or other liquids. Therefore, use the device only in dry interior rooms (e.g. never in the bath, sauna).

Never try to retrieve the device if it falls into water. Disconnect the plug immediately from the mains outlet.

Do not use the device if damage to the device itself or any of its accessories is evident.

Do not subject the device to any impact and do not drop it.

Repairs

WARNING

- Repairs to electrical appliances must only be carried out by qualified persons. Incorrect repairs can lead to considerable danger for the user. In the event of repairs, please contact our customer service or an authorized dealer.
- The zip fastener on the Shiatsu massage device must not be opened. It is only fitted for reasons relating to the manufacturing of the device.
- If the power supply unit or cord is damaged, it must be replaced by an authorised dealer.

WARNING: Fire hazard

There is a risk of fire if the device is used incorrectly or if it is not used according to these operating instructions.

For this reason, the massage device must only be operated as follows:

- Never use the device under covers, such as blankets or pillows.
- Never use the device in the vicinity of petrol or other easily flammable materials.

Handling

CAUTION

Switch the device off and disconnect the power supply unit after every use and before every cleaning.

- Never insert any objects into openings on the device, and never insert anything into the rotating parts. Make sure that all moving parts can move freely at all times.
- Only operate the device when it has been completely assembled (see chapter on Getting started).
- Never sit or stand on the moving parts of the device with all your weight and do not place objects on the device.
- Keep the device away from high temperatures.

Disposal


CAUTION

Follow local regulations regarding the disposal of materials. The appliance should be disposed of according to Regulation 2002/96/EC-WEEE (Waste Electrical and Electronic Equipment). In case of queries, contact the municipal authorities responsible for waste disposal in your area.


5. Appliance description

Shiatsu Massage Chair (front)

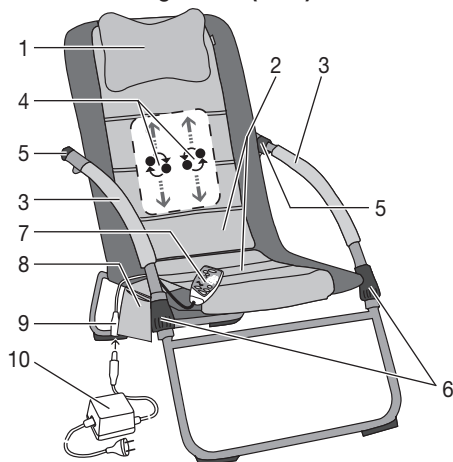- 1 Neck pillow
- 2 Removable protective covers (seat, back)
- 3 Removable protective covers (arm rests)
- 4 Four massage heads
- 5 Set the position of the back rest
- 6 Knob to adjust the seat height
- 7 Hand switch
- 8 Storage bag
- 9 Connecting coupling
- 10 Power supply unit with plug connector and power cable

Hand switch

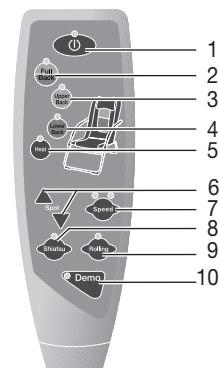- 1 Button 
- 2 FULL BACK button (entire back massage)
- 3 UPPER BACK button (upper back massage)
- 4 LOWER BACK button (lower back massage)
- 5 HEAT button (heat)
- 6 SPOT buttons (localised massage)
- 7 SPEED button (massage intensity)
- 8 SHIATSU button (Shiatsu massage)
- 9 ROLLING button (rolling massage)
- 10 DEMO button

6. Start-up

CAUTION

- **IMPORTANT:** Before using the device for the first time ensure you have released the transport lock for the massage mechanism. To do this, remove the screw using the tool provided on the rear of the massager. If the transport lock screw is not first removed and the device is used straight away, the device will become irreparably damaged and unusable. The warranty will not apply in this instance.
- Reuse the transport lock screw when transporting the device in the future. Transporting the device without using the transport lock could cause irreparable damage to the device.
- Assemble the chair by resting it on the rear legs and unfolding it.
- Back rest and seat height can be adjusted individually.
 - Seat height: To adjust the seat height, turn the knobs at the same time and make sure that the two knobs' end positions are at the same height.
 - Back rest: Hold on to the back rest and turn the ring (5) until the stopper fits through the opening and allows you to adjust the back rest. When you have reached the desired position, turn the opening of the ring away from the stopper. Complete this procedure consecutively for both sides and make sure that the end position of the rings is identical on both sides.

- **There must be no weight placed on the chair while you are adjusting seat height or back rest.**
- Secure the protective covers on the arm rests. Please make sure that you remove these before you fold the chair together again.
- Connect the power supply with the massager using the connector plug.
- Connect the device to the mains by plugging the plug into the socket.

7. Operation

During first use, place the back rest in the most upright position possible.

Gently press your back against the massage device to start with. Then make sure that you are in a comfortable position before gradually shifting your weight onto the massage device.

The massage should feel pleasant and relaxing at all times. Stop the massage or change its position, the position of the back rest or the seat height or contact pressure if you find the massage painful or unpleasant.

- Do not force any part of your body, in particular your fingers, between the rotating massage heads and their mountings in the device.
- Sit comfortably on the massage chair and make sure that you lean against it with your entire back, that you sit centrally and that the massage heads can move to the left and right of your spine.
- Switch on the massager with the  button.

The device starts automatically with the Shiatsu massage on the full back setting. The LED shows the current massage region.

You can select other settings as described below:

1. Shiatsu massage

Press the SHIATSU button. Select the desired massage region with FULL BACK (button 2), UPPER BACK (button 3) and LOWER BACK (button 4).

Spot massage, i.e. Shiatsu massage without vertical movement, can be selected by pressing one of the two SPOT buttons (number 6). The position of the localised massage can be changed by pressing and holding one of the two SPOT buttons until the desired position is reached.

Message intensity

When using the Shiatsu Massage, two different massage speeds are available, which can be selected using the SPEED button. If you have the rolling massage activated, first press the SHIATSU button in order to change to Shiatsu Massage and then set the desired massage speed.

2. Rolling massage

Press the ROLLING button. Select the desired massage region with FULL BACK (button 2), UPPER BACK (button 3) and LOWER BACK (button 4).

When the massage type (Shiatsu or rolling massage) is changed, the preselected region (full, upper, lower back) is retained. Note that the massage region must be selected again when switching from localised Shiatsu massage to rolling massage.

Heating function

Press the HEAT button (button 5) to switch on the heating function of the massage heads. The blue lights in the massage heads change colour to red.

- Use the massager for max. 15 minutes and then leave it to cool down. Long massages can lead to overstimulation of the muscles, creating tension instead of aiding relaxation. You can end the massage at any time by pressing the button . Please note that the massage heads must first return to their start position when the device is switched off. This is indicated by a flashing LED light. Once the position is reached, the device switches off.
- If the device is stopped by pressing a massage area button twice, it also returns to the starting position after 15 minutes and then switches off.

- Do not disconnect the device from the mains supply during operation. Always switch it off first and wait until it has returned to the initial position.
- Disconnect the power supply unit plug from the mains after using the device.

Note

- The device has an automatic shut-off mechanism set to a maximum operating period of 15 minutes.

8. Care and storage

CLEANING

WARNING

- Before cleaning disconnect the power supply unit from the mains socket and disconnect it from the device.
- Only clean the appliance in the manner specified. Fluid must never be allowed to enter into the device or accessories.
- Smaller stains can be removed using a cloth or moist sponge together with a little liquid detergent for delicate fabrics if necessary. Do not use any cleaning agents containing solvents.
- The removable covers can be washed at 30°C. Please note the symbols for care and cleaning on the label sewn into the cover.
Remove the filling from the neck cushion before washing it.
- Do not use the device again until it is completely dry.

Storage

If you do not use the device for longer periods, we recommend storing it in the original packaging in a dry place without weighting it down.

CAUTION

- You can fold the device for storage. Before doing so, remove the arm rests, place the back rest in the most upright position and turn the seat height to the very lowest position.
- If you use the screw for transport locking, take care when using the device again that it is definitely removed before switching on.

9. In the event of problems

Problem	Cause	Remedy
Massage heads rotating at reduced speed.	The load on the massage heads is too great.	Reduce the gradient of the back rest. Reduce the contact pressure.
Massage heads only rotate in the upper or lower region.	The massage has been activated for the limited region.	Press the FULL BACK, UPPER BACK or LOWER BACK buttons to change the massage region.
Massage heads not moving.	Appliance not connected to the power supply.	Connect the plug to the mains supply and switch on the appliance.
	Device not switched on.	Switch the device on with  . Select the massage region.
	Overheating protection tripped, device is damaged.	Contact your authorized dealer or the service address.

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