beurer MG 170

- D** **Massage-Matte**
Gebrauchsanleitung
- GB** **Massage mat**
Instruction for Use
- F** **Tapis de massage**
Mode d'emploi
- E** **Colchón de masajes**
Instrucciones para el uso

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Items included in the package:

- Massage mat
- Power supply unit
- These operating instructions

Explanation of symbols

The following symbols are used in the operating instructions and on the type plate:

WARNING Warning of risks of injury or health hazards.

CAUTION Safety information about possible damage to appliance/accessories.

Note Important information.

The device has double protective insulation and therefore complies with protection class 2.

Only use indoors.

Please read the instructions!

Do not use folded or rolled up!

Do not pierce with needles!

1. Getting to know your appliance

Instrument functions

You can massage your entire body with this electric massager. Massages can be both relaxing and stimulating, and are often used in the treatment of stiff or tense muscles, pains and fatigue. This device massages with the aid of 8 vibration motors and offers a wide range of setting options.

The device is only for **private use** and is not intended for medical or commercial purposes.

2. Proper use

This device is exclusively designed for massaging the human body (e.g. back, legs). It cannot replace medical treatment. Do not use the massage device if any of the following warnings apply to you.

WARNING

Use of the massage device is prohibited in the following cases:

- Not when driving a motor vehicle,
- not with changes caused by disease or injury in the back region (e.g. disc problems, open wounds),
- not on children,
- not during pregnancy,
- never while you sleep,
- not on animals,
- not during activities during which an unexpected reaction may be dangerous,
- not while ingesting substances that reduce reaction times (e.g. painkiller medications, alcohol),
- never longer than 15 minutes (risk of overheating) and allow it to cool off for at least 15 minutes before using it again.

Please consult your doctor if you are unsure whether or not the massage device is suitable for you.

This device is not intended for use by persons (including children) with restricted physical, sensory or mental capabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety or they are given instructions on how to use the device.

Children should be supervised to ensure that they do not use the device as a toy.

The device has a hot surface in the area of the heating element. Persons who are sensitive to heat must be careful when using the device.

You should consult your doctor before using the massage device if:

- if you are suffering from a serious illness or have undergone an operation in the area to be massaged,
- you are fitted with a cardiac pacemaker, implant or other aid,
- Thromboses,
- Diabetes,
- in case of pains of unclarified origin.

The device is only intended for the purpose described in these operating instructions. The manufacturer cannot be held liable for damage occurring as a result of inappropriate or careless use.

3. Safety information

Please read these operating instructions carefully. Failure to observe the information provided below can lead to personal injury or material damage. Keep these instructions in a safe place where they can also be read by other persons using the device. Always pass on these instructions together with the appliance if you sell it or pass it on.

WARNING

Keep children away from the packaging materials. Risk of suffocation.

WARNING: Electric shock

Just like any other electrical appliance, this massage device must be handled with due care and attention in order to prevent the risk of electric shock.

For this reason, the device must only be operated as follows:

- Only with the supplied power supply unit at the mains voltage indicated on the power supply unit.
- Do not use the device if damage to the device itself or any of its accessories is evident.
- Never use the appliance during a thunderstorm.

In the event of a defect or malfunction, switch the device off immediately and disconnect it from the power supply. Never pull on the power cable or on the appliance in order to disconnect the power supply from the mains outlet. Never hold or carry the device by the power cable. Keep the cables away from hot surfaces. Do not pinch, kink or twist the cables. Do not stick in any pins, needles or sharp objects.

Make sure that the massage device, the switch, the power supply unit and the cable do not come into contact with water or other fluids. For this reason, the device must only be operated as follows:

- Only in dry inside rooms.

Never try to retrieve the device if it falls into water. Disconnect the plug immediately from the mains outlet.

Do not use the device if damage to the device itself or any of its accessories is evident.

Do not subject the device to any impact and do not drop it.

Repairs

WARNING

- Repairs to electrical appliances must only be carried out by qualified persons. Incorrect repairs can lead to considerable danger for the user. In the event of repairs, please contact our customer service or an authorized dealer.
- The zip fastening of the massager must not be opened. It is only used for technical production reasons.
- If the power supply unit or its lead is damaged, they must be replaced by an authorised retailer.

Fire hazard

WARNING

There is a risk of fire if the device is used incorrectly or if it is not used according to these operating instructions.

For this reason, the massage device must only be operated as follows:

- Never use the device under covers, such as blankets or pillows.
- Never use the device in the vicinity of petrol or other easily flammable materials.

Handling

CAUTION

Switch the device off and disconnect the power supply unit after every use and before every cleaning.

- Keep the device away from high temperatures.

Disposal

CAUTION

Follow local regulations regarding the disposal of materials. The appliance should be disposed of according to Regulation 2002/96/EC-WEEE (Waste Electrical and Electronic Equipment). In case of queries, contact the municipal authorities responsible for waste disposal in your area.

4. Start-up

- Remove all packaging and transport locks.
Ensure that the packaging foils do not fall into children's hands. Risk of suffocation.
- Check the device, power supply unit and all cables for damage.
- Position the massage mat on a suitable surface (e.g. floor, bed, couch).
- In order to secure the mat, use the flexible tie straps if necessary.
- Connect the power supply unit to the mat using the connector plug.


WARNING

- Do not subject the device to any impact and do not drop it.
- Do not use the device if damage is evident to the device itself or any of its accessories.
- Use the device only with the power supply unit provided and using only the voltage specified on the power supply unit.


5. Appliance description

Massage mat (front)

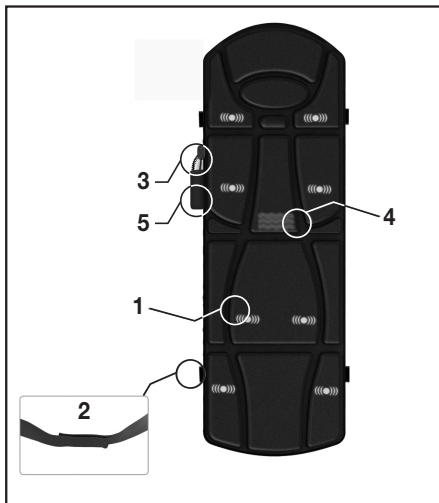- 1 8 massage motors
- 2 Fastening straps
- 3 Hand switch
- 4 Heat
- 5 Storage bag

Hand switch

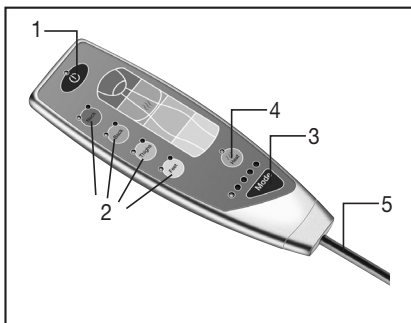- 1 ① button with LED indicator
- 2 Massage areas 1-4 (can be individually selected), each with 2 intensity levels
- 3 MODE button with LED indicators
- 4 🌊 button with LED indicator
- 5 Connection cable to the device

6. Operation

- Lie or sit down comfortably on the mat.
- Switch on the massage device with the ① button.
- Using MODE, choose the desired massage program.
- Using buttons 1–4, select the area to be massaged and the desired massage intensity. The activated massage zones will be shown by the LED indicators.
- The automatic switch-off function is activated by starting the massage and is preset to 15 minutes.
- Press the 🌊 button to switch the temperature function on and off.
- Do not use the massage device for longer than 15 minutes. The massage can be finished at any time by pressing the ① button.
- After the massage is finished, disconnect the device from the mains.

**Note**

- Do not use the massage device just before going to bed. The massage can also have a stimulating effect and may cause difficulties in getting off to sleep.

**WARNING**

The massage must be pleasant and relaxing at all times. Stop the massage or change your position or the contact pressure if you find it is painful or unpleasant.

7. Care and storage

CLEANING

**WARNING**

- Disconnect the device before cleaning it.
- Only clean the appliance in the manner specified. Fluid must never be allowed to enter into the device or accessories.
- Do not use the device again until it is completely dry.

First unplug the device and allow it to cool off. Then you can clean it with a cloth which is no more than lightly moistened. You can also use a mild detergent.

Storage

- Keep the massage device in a cool and dry place where it is out of the reach of children. We recommend storing the device in its original packaging.
- Avoid any contact with sharp or pointed objects which could damage the device or cover.

8. In the event of problems

Problem	Cause	Remedy
The device is very loud during massage.	The seat cover is not fully in The mat is not in full contact with the surface it was placed on.	Ensure that the entire mat is resting flat
Massage stopped.	The preset maximum period of operation has elapsed.	Let the device cool down for at least 15 min. before using again. Switch the device on again.
	The connection to the mains has been interrupted.	Check the plug connection and the contact of the mains plug with the mains supply.
No massage.	Appliance not connected to the mains power supply.	Make all connections and switch on the device.
	Device not switched on.	Turn the device on using  and activate the massage using buttons 1–4.

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