beurer MG 140

- (D)** Shiatsu-Massagekissen
Gebrauchsanleitung
- (GB)** Shiatsu massage cushion
Instruction for Use
- (F)** Coussin de massage Shiatsu
Mode d'emploi
- (E)** Cojín de masaje Shiatsu
Instrucciones para el uso
- (I)** Cuscino massaggiante shiatsu
Istruzioni per l'uso
- (TR)** Shiatsu masaj yastığı
Kullanma Talimatı
- (RUS)** Массажная подушка шиatsu
Инструкция по применению
- (PL)** Shiatsu – poduszka do masażu
Instrukcja obsługi

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Items included in the package:

- Shiatsu massage pillow
- Removable cover
- Power supply unit
- These operating instructions

Explanation of symbols

The following symbols are used in the operating instructions and on the type plate:

WARNING Warning of risks of injury or health hazards

CAUTION Safety information about possible damage to appliance/accessories

Note Important information

The device has double protective insulation and therefore complies with protection class 2

1. Getting to know your appliance

Instrument functions

With the aid of this electrical Shiatsu massage device, you can give yourself or somebody else an effective and pleasant massage. Massages can be both relaxing and stimulating, and are often used in the treatment of stiff or tense muscles, pains and fatigue. The device has a universal wedge shape which enables it to deliver powerful and intensive deep massages for neck, back, legs and arms. The device can also be placed on the floor for an invigorating foot massage.

Thanks to its compact size and low weight, the device is easy to transport and put away when not in use.

The device is for **private use only and is not intended** for medical or commercial purposes.

In addition, the Shiatsu massage device offers the following advantages:

- Choice of two massaging directions
- Simple operation
- Covers can be removed for washing
- Power supply unit can be disconnected

2. Proper use

This device is intended only for the massaging of individual parts of the human body (see chapter 6). It cannot replace medical treatment. Do not use the massage device if any of the following warnings apply to you. Please consult your doctor if you are unsure whether or not the massage device is suitable for you.

WARNING

Use of the massage device is prohibited in the following cases:

- On children under the age of 15
- On animals
- On persons suffering from a herniated disk or from abnormal changes to the vertebrae
- Never use the massage device to massage the heart region of persons fitted with a pacemaker. Always consult a doctor before using the massage device on other parts of the body of a person fitted with a pacemaker.

- Never use the massage device on swollen, burned, inflamed or damaged skin or body parts.
- The massage device must not be used on lacerated wounds, capillaries, varicose veins, acne, couperose, herpes or other skin diseases.
- The massage device must not be used during pregnancy.
- Do not use the massage device on the face (eyes), on the larynx or on other particularly sensitive parts of the body.
- Never use the device while you are asleep.
- Never use the device after taking medication or consuming alcohol (reduced levels of awareness).
- Do not use the device in the car.

This device is not intended for use by persons (including children) with restricted physical, sensory or mental capabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety or they are given instructions on how to use the device.

Children should be supervised to ensure that they do not use the device as a toy.

You should consult your doctor before using the massage device if:

- you suffer from a serious illness or are recovering from an operation;
- you suffer from any of the following: common cold accompanied by a fever, varicose veins, thrombosis, phlebitis, jaundice, diabetes, nervous system diseases (e.g. sciatica) or acute inflammations.

The device is only intended for the purpose described in these operating instructions. The manufacturer cannot be held liable for damage occurring as a result of inappropriate or careless use.

3. Safety information

Please read these operating instructions carefully. Failure to observe the information provided below can lead to personal injury or material damage. Keep these instructions in a safe place where they can also be read by other persons using the device. Always pass on these instructions together with the appliance if you sell it or pass it on.

WARNING

Keep children away from the packaging materials. Risk of suffocation.

WARNING: Electric shock

Just like any other electrical appliance, this massage device must be handled with due care and attention in order to prevent the risk of electric shock.

For this reason, the device must only be operated as follows:

- The massage device must only be operated with the supplied power supply unit, which must be connected to a power supply with the mains voltage indicated on the power supply unit.
- Never use the device if any damage to the appliance itself or to any of its accessories is evident.
- Never use the device in rooms in which aerosols are sprayed or oxygen is administered.
- Never use the appliance during a thunderstorm.

In the event of a defect or malfunction, switch the device off immediately and disconnect it from the power supply. Never pull on the power cable or on the appliance in order to disconnect the power supply from the mains outlet. Never hold or carry the device by the power cable. Keep the cables away from hot surfaces.

Make sure that the massage device, the power supply unit and the cable do not come into contact with water or other fluids. For this reason, the device must only be operated as follows:

- Only in dry inside rooms
- Only with dry hands
- Never in the bath, in the shower or near a filled basin
- Never in swimming pools, whirlpools or in the sauna
- Never out in the open

Never try to retrieve the device if it falls into water. Disconnect the plug immediately from the mains outlet.

Repairs

WARNING

- Repairs to electrical appliances must only be carried out by qualified persons. Incorrect repairs can lead to considerable danger for the user. In the event of repairs, please contact our customer service or an authorised dealer.
- The zip fastener on the Shiatsu massage device must not be opened. It is only fitted for reasons relating to the manufacturing of the device.

WARNING: Fire hazard

There is a risk of fire if the device is used incorrectly or if it is not used according to these operating instructions.

For this reason, the massage device must only be operated as follows:

- Never operate the massage device for longer than 15 minutes (risk of overheating) and allow it to cool off for 15 minutes before using it again.
- Never leave the device unsupervised, particularly if children are around.
- Never use the device under covers, such as blankets or pillows.
- Never use the device in the vicinity of petrol or other easily flammable materials.

Handling

CAUTION

Switch the device off and disconnect the power supply after every use and before every cleaning.

- Never insert any objects into openings on the device, and never insert anything into the rotating parts. Make sure that all moving parts can move freely at all times.
- Make sure that no hair can get caught in the device. Wear a hair net if necessary.
- Avoid contact with sharp or pointed objects.
- Do not kink the device or cable.
- Never sit or stand on the device with all your weight and do not place objects on the device.
- Do not place the device on delicate surfaces, as otherwise the zip fastener could potentially damage the surface.
- Make sure that the appliance is not too close to heat sources such as ovens, radiators or direct sunlight.

Disposal


CAUTION

Follow local regulations regarding the disposal of materials. The appliance should be disposed of according to Regulation 2002/96/EC-WEEE (Waste Electrical and Electronic Equipment). In case of queries, contact the municipal authorities responsible for waste disposal in your area.

4. Appliance description

Shiatsu massage pillow (front)

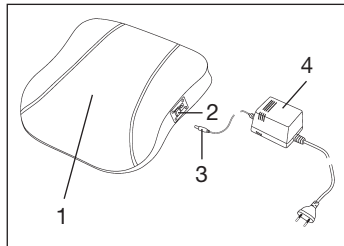- 1: Four massage heads, rotating in pairs
- 2: Switch, 3-level
 - Level 1: Switch on UPWARD massage
 - Level 0: Switch message device OFF
 - Level II: Switch on DOWNWARD massage
- 3: Plug inserted in socket (can be pulled out)
- 4: Power supply unit

5. Start-up

- Remove all packaging and transport locks.
- Check the device, power supply and cables for damage.
- Fit the protective cover to the appliance if it is not already in place.
- Connect the power supply plug to the massage device. The switch must be set to "0". Move it to "0" if it is not already in this position.

WARNING

- Do not use the device if damage is evident to the device itself or any of its accessories.
 - The massage device must only be operated with the supplied power supply unit, which must be connected to a power supply with the mains voltage indicated on the power supply unit.
 - Route the cable safely so that there is no risk of anyone tripping over it.
- Connect your device to the mains by plugging the power supply unit into the mains outlet.

6. Operation

Use the switch on the side to switch on the massage device. The switch positions „I“ and „II“ also determine the massaging direction:

- Switch position I: the massage heads rotate against each other in an upward motion.
- Switch position II: the massage heads rotate against each other in a downward motion.

Move the switch to the middle position (level 0) to switch the device off.

Note

- Always wear light clothing when receiving a massage from the device.
- To receive a more gentle massage, insert a folded towel between the massage device and your body.
- Do not use the massage device for longer than 15 minutes. Otherwise, using the device for longer can over-stimulate the muscles and cause tension instead of relaxation.
- Do not use the massage device just before going to bed. The massage can also have a stimulating effect and may cause difficulties in getting off to sleep.

WARNING

The massage must be pleasant and relaxing at all times. Stop the massage or change the position of the massage device if you find it is painful or unpleasant.

- In order to prevent yourself from being trapped between the heads, do NOT lean on the massage heads when changing your position or moving the massage device.
- Do not force any part of your body, in particular your fingers, between the rotating massage heads and their mountings in the device.


To start with, gently press the part of your body you wish to massage against the massage device. Particularly when lying down, make sure that you gently lift your weight off the corresponding part of your body. Then make sure that you are in a comfortable position before gradually shifting your weight onto the massage device.


Back massage: upper back area


Back massage: lower back area (lumbar vertebrae)

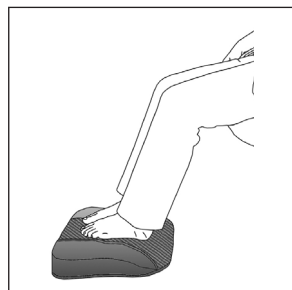Foot massage

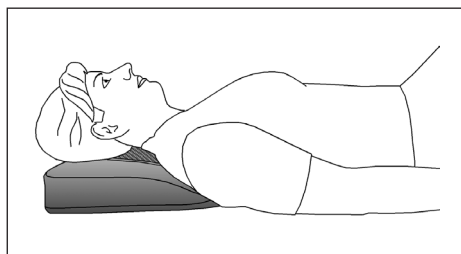Neck massage

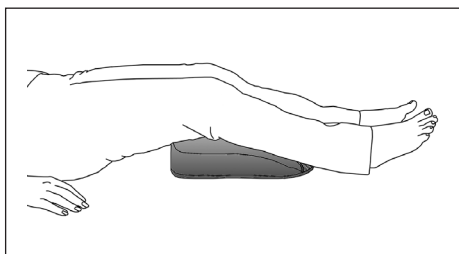Leg massage

7. Care and storage

Cleaning

WARNING

- Disconnect the device before cleaning it.
- Only clean the appliance in the manner specified. Fluid must never be allowed to enter into the device or accessories.
- Do not use the device again until it is completely dry.

First unplug the device and allow it to cool off. Then you can clean it with a cloth which is no more than lightly moistened. You can also use a mild detergent. The removable cover is machine-washable at 40 °C. Please note the cleaning and care symbols on the label which is sewn into the cover.

Storage

- Keep the massage device in a cool and dry place where it is out of the reach of children. We recommend storing the device in its original packaging.
- Avoid any contact with sharp or pointed objects which could damage the device or cover.

8. In the event of problems

Problem	Cause	Remedy
Massage heads rotating at reduced speed.	The load on the massage heads is too great.	Release the massage heads.
Massage heads not moving.	Appliance not connected to the power supply.	Connect the plug to the mains supply and switch on the appliance.
	Switch set to "0".	Switch on the appliance (switch position I or II).
	Overheating protection tripped.	Disconnect the device from the power supply, wait at least 15, release any blocked massage heads, reconnect and switch on.

9. About Shiatsu massage

Shiatsu is a form of body massage which was originally developed in Japan around 100 years ago. It is based on ideas anchored in traditional Chinese medicine, in particular the meridian system, i.e. the paths in the human body along which energy is conveyed.

The aim of any Shiatsu massage is to promote the physical, emotional and mental wellbeing of the person being treated. To do this, energy blocks and bottlenecks in the meridians need to be dissolved and the self-regulating ability of the body needs to be stimulated.

Shiatsu therapists achieve this by exercising pressure in flowing movements along the meridians. Although Shiatsu can be literally translated as "finger pressure", practitioners also use the heels of their hands, their elbows and knees as well as their fingers. In the process, the therapist will use his or her full body weight in order to generate the required pressure. Your Shiatsu massage device has rotating massage heads which imitate the pressure and kneading motions of a Shiatsu massage.

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