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Give nature a helping hand!

By purchasing the Original **NOVAFON®** instrument you have definitely taken a positive step in the right direction for your health. These intrasonic wave instruments are specially designed for areal and selected application.

List of Contents

1. Technical Specifications
2. Instrument Description
3. Operating Instructions
4. Operating Mode
5. Fields of Application
6. Examples of Use
7. Contraindications
8. Using the magnetic probe.
Sonic and magnetic therapy
9. Sole of the Foot – Reflex Zones – Massage
10. Warranty
11. Customer Service

1. Technical Specifications

- ~ Alternating Current
- OFF (a part of the instrument only)
- ◉ ON (a part of the instrument only)
- A Protection Class II Device
- hPA** Hecto Pascal
- Hz** Hertz
- KB** Short-Time Operation
- mA** Milliampere
- V** Volts

Product Information:

Category of Product: Electro-medical
massage and nerve stimulation device.

Mains supply voltage: primary 230 V~ 50 Hz 65 mA,
secondary 17~ 50 Hz 450 mA

Nominal Frequency: 50 Hz

Nominal Intake: 9.5 W

Power Supply: STO 15–20 mins

Protection Classification: II

Manufacturer: NOVAFON® GmbH

Daimlerstraße 13, D-71384 Weinstadt, Germany

2. Instrument Description

SK1/1 chromed

This particular model is the original **NOVAFON®** sonic wave instrument in its earliest form. It operates using sound waves or vibrations produced by applying an electrical current to magnetic trip coils. The resulting sonic waves are applied vertically to the body. The SK1/1 model comes equipped with a single switch setting producing 120 vibrations per second.

SK2/2 chromed

This model is a further development of the type SK1/1. This model comes equipped with a dual switch setting. It includes:

- Setting 1 = normal mode
- Setting 2 = percussion mode

Setting 1 : 100 – 120 vibrations per second

Setting 2: 50 – 60 vibrations per second

Using Setting 2 the operator is able to perform with ease a precisely applied pummelling massage and acupuncture treatment. This drumming effect can be used each time treatment is given provided there is no discomfort involved. The retractable, extendable handle ensures that even the most difficult areas of the body, like shoulders, back, back of the neck etc., can be reached with ease.

SK1/1 plastic

Made of plastic, the SK1/1 model was devised for cosmetic purposes. It too works using sound waves, however, with a reduced penetrative effect. Similar to the chromed-plated SK1/1 model, this version has only one setting, vibrating 120 times per second and is not equipped with an extendable handle. Each model comes complete with two interchangeable fittings: a sonic plate and a sonic ball. Your original **NOVAFON®** instrument is fitted with a thermal protection device which, in case of overheating, switches off the appliance thus preventing damage.

3. Operating Instructions

Initial Operation

Insert the plug into a power point with the corresponding electrical current as specified above. Ensure that the power cable is not stretched beyond the advised 1.5 metres otherwise there is a risk of breakage. Use an extension lead if necessary.

SK1 model

Turn off: OFF ○

Turn on: ON ⊙

SK 2 model (with percussion mode)

middle setting: OFF

- 1 = standard operating mode
- 2 = percussion mode

Power control

To regulate the power turn the rotatable control dial. Viewed from above, rotating the dial to the right in a clockwise direction decreases the intensity. Correspondingly, turning the dial to the left in an anti-clockwise direction increases the intensity. By adjusting it in this way, you can set the appliance to your own personal needs. It is recommended beginning with a low amplitude, then increasing it if and when necessary. Avoid acting according to the adage 'more helps more'. That would be misguided. Select the intensity with which you feel comfortable. Then your muscles will feel the same. Applying excessive pressure to the appliance does not increase its effectiveness, but instead impedes its vibratory performance. The optimum vibration is achieved solely by means of the appliance's own bearing-strength. Sonic plate and sonic ball Both attachments can be interchanged by simply unscrewing whilst applying light pressure to the device. In accordance with EN 30993, these parts contain no cytotoxicants. Should these parts need replacing due to normal wear, then you can order them from your dealer accordingly. Further important information concerning use The original **NOVAFON®** instrument is an electro-medical instrument and must be handled and stored with care.

- Store only in a dry, clean place at a temperature between + 10° C and + 40° C
- Connect only to the correct alternating current power supply and ensure that your mains voltage corresponds to the voltage specifications stated on the transformer.
- Use a suitable extension cable when further than 1.5 metres away from the power point.
- Do not bring the instrument into contact with water or other liquids.
- Use both the appliance and the transformer for no more than 20 minutes at any given time, allowing at least a 15 minute break in order for them to cool down sufficiently.
- Avoid banging or dropping the appliance.
- Only use creams, salves or other cosmetic and medical ointments according to instructions. Deposits of such compounds may damage the instrument. We recommend the use of a fine cloth placed between the instrument and the skin if creams or ointments have been applied prior to use. Using the instrument before applying such compounds helps guarantee the circulatory enhancing effect required by tissue to absorb the substances they contain.
- For hygienic reasons, make sure that the respective utensils are not used on more than one person at any one time.
- If you have purchased an additional magnetic probe, then it should be kept away from delicate electronic equipment (e.g. clocks and watches) and magnetic recording devices such as credit and debit cards as the strong magnetic field may partly or totally erase any data recorded.
- During treatment you should hold the Original **NOVAFON®** instrument by the handle or retractable extension only.
- In order to clean the instrument's outer casing and additional fittings we recommend using a clean, damp cloth.

Damaged transformers and equipment should be withdrawn from use.

The NOVAFON® instrument is an electromagnetic device which must only be used in accordance with the purposes outlined in this leaflet or its related publications.

4. Operating Mode

For many years now doctors, therapists and clinics have successfully been making use of the Original **NOVAFON**[®] instrument in their treatment of various forms of medical complaints. Used by patients themselves, it has proved itself to be a positive, supporting and continuing treatment adopted in the remedy of a variety of chronic and mild complaints as well as in preventative care.

Sound waves

The **NOVAFON**[®] instrument produces intrasonic waves in a mix frequency of between 100 and 400 Hertz (main frequency 100 Hz). Thus, the frequency span is within audible range, the frequency being the amount of vibrations per second. Healthy human cells also vibrate within audible range. As a result, a feature of the **NOVAFON**[®] instrument is to actively stimulate the cells' own natural vibration, and that up to a depth of 60 mm or more.

Cells

The active stimulation of the cells' own natural vibration along with the instrument's mechanical pounding effect have an ideal impact upon cellular metabolism. Consequently, intercellular activity is improved due to the increased permeability of the cell membrane; the diffusion process is regulated; and substances otherwise harmful to cells are dispersed more quickly. Due to the tissue's reactively stimulated circulation, the molecules in that tissue are nourished with increased levels of oxygen and nutrients.

Nerves and muscles

Moreover, vibration can lower the perception threshold to touch, pain and warmth, thus explaining the pain-soothing effect of the **NOVAFON**[®] instrument (stimulation of the nerve endings). Above all, it is the regulating effect upon muscle tension – whether too high or too low – due to stimulation of the receptors of both muscles and tendons which is the reason why the **NOVAFON**[®] instrument can be used with such versatility.

Sonic plate and sonic ball

Whenever you treat a flat, easily massaged part of the body, apply the sonic plate. It should lie in a flat position on the skin. For parts of the body that are less easily reached – such as between the fingers, toes and ribs, the area around the armpits or behind the ears – the sonic plate is best used. The ball is also utilised for acupressure and to massage reflex zones.

5. Fields of Application

For a long time now the **NOVAFON**[®] instrument has been used successfully to treat ailments of all kinds. It is passed over the area to be treated using either linear strokes or slowly rotating movements. Where possible, the treatments should be performed at least twice a day and should last no longer than 15 to 20 minutes.


General application ranges are:

Shoulders

Neck and Back

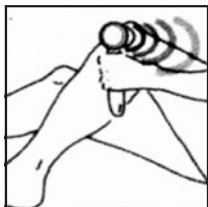
Knees, lower leg and ankles

Lower back and Abdomen

Elbows and Arms

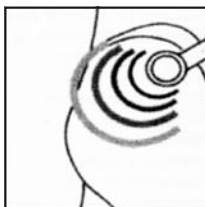Hips, buttocks and thighss

When treating chronic ailments, it should be assumed that these will initially revert back to an acute state so that the body's natural regenerative capability can be activated. This is referred to as an initial deterioration. Then, experts also say, that once the body's regenerative capability has been awakened, continuing with the treatment leads to success. Pain and everyday ailments can also be signs warning of serious illness. It is absolutely necessary that the cause be treated and it is therefore essential that you undergo professional medical care. Making use of the **NOVAFON®** instrument is however a worthwhile aid that can also assist your doctor in prescribing complementary, physical therapeutic measures. Such therapeutic measures are prescribed where there is pain and complaints of the joints (incl. muscles and tendons), sport injuries and their prevention, including:

Achilles tendon complaints, Arthrosis, Bechterew's disease, Connective tissue disorders, Circulatory and sensory disorders, Tennis elbow, Fractures (aftercare), Joint/muscle fatigue, Arthropathy, Lumbago, Pain in the hip joint, Cervical spine syndrome, Sciatica, Knee injuries, Bone injuries, Lower back pain, General massage, Muscular disorders, Myogelosis, Neuritis, Trapped nerve, Neurological disorders, Paresthesia, Contusions and hematoma, Rheumatism of the joints, Dislocated shoulder, Tendonitis, Sport massage, Ganglion, Sprains, Dislocations, Leg cramp, Spinal blockages, Flail and instable joints, Strained muscles

When using the **NOVAFON®** instrument – as always when undertaking any physical treatment – the most important rule is: always avoid irritable joints. Treat around them. The following can also be treated without constant medical supervision at home:

Effects of ageing, Bumps and protrusions, Bee stings (insect bites), Chronic pain (like migraine), Constipation, Colds, Wrinkles and lines, Heel spur Foot treatment, Bronchitis, Circulation problems, Lipoma (benign), Lymphostasis, Menstrual problems, Keloid scars, Sinus problems & sinusitis, Edema (post operative), Cellulite, Vegetative imbalance, Generalized tendomyopathy

Further therapeutic uses and more information about possible forms of treatment along with their effects are to be found in our selected range of publications (see page 24) as well as on our homepage www.novafon.de where reports by doctors and therapists pertaining to the use of the **NOVAFON®** instrument are readily available.

6. Examples of Use

Headache

Headaches can have many different reasons. The root cause of persisting headaches should be diagnosed by a physician. Often it is a case of tension in the muscles of back of the neck leading to a lack of blood circulating to the head thus triggering headache. In such cases apply the NOVAFON® instrument several times a day for a maximum of 5 minutes as shown in Figure 1.

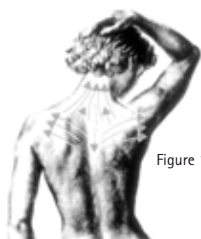
Figure 1

Arthrosis of the hip

Include in the treatment all muscles surrounding the hip joint, that is, those belonging to the thigh, buttocks and the base of the spine as depicted in Figure 2. As the unaffected side often compensates for the weight of the infected joint, it may well be overburdened. It therefore makes sense to treat the healthy side too.

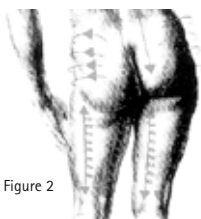
Figure 2

Tenosynovitis (Inflammation of the tendons)

Tendons sheaths consist of fibrous tissue in which the tendons glide back and forth. In the case of inflammation, treat the entire forearm from the wrist up to the elbow in a circular motion, and this at best mornings and evenings (Fig. 3)


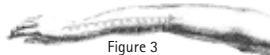
Figure 3


Tennis elbow

So-called tennis elbow is an inflammation of the elbow joint. In many instances, this chronic sprain is accompanied by an inflammation of the articular capsule and/or of the synovial bursa, bone membrane (periosteum), tendons and muscles. Treat as shown in Figures 4, 4a and 4b, keeping mainly to the muscular tissue and avoiding close proximity to the bone, where it is otherwise possible to aggravate the bone membrane.

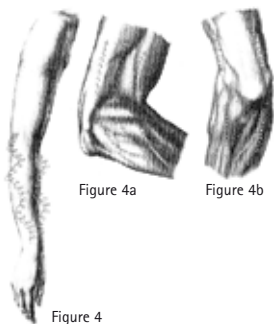
Lumbago and Sciatica

Distinctive symptoms of lumbago are sudden, shooting pains in the lower back or lumbar (for treatment see Fig. 5). Different from lumbago, the symptoms of sciatica mainly include a contracting pain radiating from the lumbar region over the buttocks into the left or right leg. Invariably, acute sciatica belongs in the hands of a physician, then there are types which make a quick operation unavoidable. Once the acute inflammation has subsided, the **NOVAFON**® instrument can be used. Treat the lumbar region of the spine along with the affected leg using rapid, circular movements.

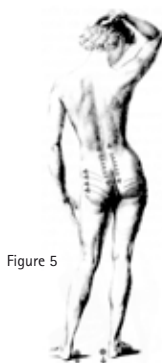
Knee injuries

When treating strained or overstretched ligaments connecting and supporting the knee joint on the inner leg, proceed with the **NOVAFON**® instrument adopting either stroking or circular movements to and fro along the lines depicted in Fig. 6.

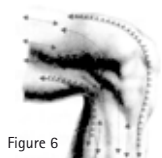Bronchitis

With bronchitis especially, the **NOVAFON**® instrument has a dual benefit. Firstly, those respiratory muscles that are normally subject to cramp are relaxed. Secondly, the vibrations entering the lungs via the ribcage act like a thumping massage thus freeing congested phlegm. Work the **NOVAFON**® instrument 'step fashion' in arch-shaped movements from the spine towards the sternum and back again. Follow the corresponding intercostal space whilst doing so. The treatment should last at least 5 minutes and be targeted towards the amount of phlegm expectorated. In the case of deep-seated phlegm, the treatment is to be repeated after 3 to 4 hours. Treatment should not exceed 15 minutes in order to prevent overexertion.

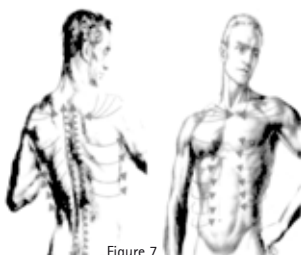7. Contraindications

The NOVAFON® instrument should definitely not be used in the following cases:

People fitted with heart pacemakers, artificial heart valves or suffering from cardiac palpitations as well as epileptics and pregnant women should consult a physician prior to use. Always avoid applying the instrument directly on or over surgical implants or inflamed joints and veins. In such cases you should also refrain from using the magnetic probe. The instrument should not be used in conjunction with inflamed diseases of the skin, dyshidrotic eczema as well as cases of inflammation often accompanied by fever or a high temperature. Also, treatment in proximity to ailing veins, and in the region of the abdomen and back of pregnant women is not recommended. At the earliest, treatment to whiplash should not commence until 3 days after the initial injury. The magnetic probe should not be applied to areas surrounding surgically-pinned fractures nor on open wounds. It is advised to wait until the fourth day before treating bruises/contusions.

If in doubt, always consult a physician.

8. Using the magnetic probe

Sonic wave therapy plus magnetic field therapy.

During magnetic field therapy, with the flow of electrically conductive blood plasma passing through the blood vessels, a vibrating (oscillating) separation of ions takes place in the area around the alternating magnetic pole fields. In addition, Faraday's law of induction comes into effect, by which electrical eddy currents occur, then subsequently release heat. Through the amalgamation of sonic wave and magnetic field therapies, the effects of both are combined. With the **NOVAFON®** instrument alone, a micro-massage of the joints' inner surfaces, in the knees, hip joints and in the joints of the arm takes place, for example. The result is a more intensive effect and an acceleration of the healing process. Of great importance to the success of this combined therapy is not only the proven depth effectiveness of the **NOVAFON®** instrument but also the gauss force of the magnets employed during that therapy. For treatment using the magnetic probe, the sonic disc is unscrewed and the magnetic probe screwed on. Coated with 24 carat gold, it is particularly friendly to the skin. The sonic therapy can also be carried out through a piece of thin material. For magnetic field therapy, skin contact is indeed essential. The relatively hard magnet should be applied gently to the skin however. All types of use, both the intrasonic therapy and the magnetic field therapy, have no known side effects when used on the body for periods of up to 20 minutes.

9. Sole of the Foot – Reflex Zones – Massage

The **NOVAFON**[®] instrument is also ideally suited to foot massage. Due to the depth reached by the sonic waves, the massage effect is more powerful and longer lasting than one achieved using conventional massage equipment or done by hand. The massage is possible using the sonic plate as well as the sonic ball. During treatment using the sonic plate, all the reflex zones are covered by employing the apparatus not only to the soles of the feet, but also to the sides and bridge of the foot. Owing to the specific effect of the intra-waves, a massage lasting approximately 2 to 4 minutes per foot achieves the optimum effect. One feels rejuvenated just like after a refreshing shower. Treatment of individual zones is done using the sonic ball or the magnetic probe. With these attachments, specifically chosen zones can be worked on, for example, those connected with parts of the body such as the shoulder, or organs like the liver. Due to the intensity of vibration, a massage carried out using the **NOVAFON**[®] instrument is not only more effective but also much less strenuous than one 'done by hand'.

10. Warranty

The appliance comes complete with a 2 year guarantee which commences on the date of transfer. The proof of purchase (i.e. receipt, bill, letter of consignment, delivery documentation etc.,) exhibits this transfer date. This should be kept in a safe place. Should you have any cause for complaint, please contact your supplier. Initially, any right of exchange, reduction or subsequent supply is excluded. You do however reserve the statutory right to claim a reduction (discount in price) or a cancellation of sale (withdrawal from purchase) should we establish a defect covered under the warranty which cannot be rectified within a reasonable period of time. Any attempt to tamper with the instrument will render the warranty null and void. Furthermore, great care must be taken when using the sound wave therapeutic instrument not to extend the power cable beyond the maximum length of 1.5 metres. Anything over and above this prescribed length requires an extension lead, as otherwise there is the risk of damage thus rendering the warranty invalid.

11. Customer Service

As the **NOVAFON**[®] instrument is a finely adjusted and precise instrument, it may only be repaired by a trained member of the **NOVAFON**[®] customer service department. However, in general, repair is more costly than buying a new instrument. Therefore, please do not return any faulty appliances to us unsolicited, but rather discuss the defect and the cost of repair with us first. If, while still under guarantee, repair is no longer cost effective, the **NOVAFON**[®] Company will offer you a replacement (once again with a 24 month warranty). Please do not dispose of the instrument in your domestic waste. The **NOVAFON**[®] Company will gladly accept the return of obsolete instruments.

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- Program 1 = normal
- Program 2 = percussion

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NOVAFON® Sound Wave Instruments GmbH