<p>(D) Massagesessel Gebrauchsanleitung</p>	2-10	<p>(I) Poltrona per massaggi Istruzioni per l'uso</p>	38-46
<p>(GB) Massage chair Instructions for use</p>	11-19	<p>(TR) Masaj koltuğu Kullanım kılavuzu</p>	47-55
<p>(F) Fauteuil de massage Mode d'emploi</p>	20-28	<p>(RUS) Массажное кресло Инструкция по применению</p>	56-64
<p>(E) Sillón de masaje Instrucciones de uso</p>	29-37	<p>(PL) Fotel do masażu Instrukcja obsługi</p>	65-73

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Dear customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of beauty, heat, gentle therapy, blood pressure/diagnosis, weight, massage and air.

Please read these instructions carefully and observe the information they contain.

With kind regards,
Your Beurer Team

1. Intended use.....	11
2. Safety notes.....	11
3. Getting to know your device	12
4. Unit description	13
5. Initial use.....	15
6. Operation.....	17
7. Maintenance and storage	18
8. What if there are problems?.....	19
9. Technical Data.....	19
10. Disposal.....	19
11. Warranty and service	19

Signs and symbols

The following symbols appear in the instructions for use and on the type plate:

Warning!

Warning instruction indicating a risk of injury or damage to your health.

Caution!

Safety note indicating possible damage to the device/accessory.

Important!

Note indicating important information.

Only use in closed rooms.

Do not insert needles!

Cross-reference to another section.

Included in delivery

- Massage chair
- Seat cushion
- Back cushion
- Pillow

1. Intended use

The massage chair is intended for the purpose of massaging the back and legs. Please observe the safety notes in these instructions for use.

Use the device only according to the methods specified. The device is suitable only for private use and not for medical or commercial purposes.

2. Safety notes

- Read these instructions for use carefully.
- Non-observance of the following information may result in personal injury or material damage.
- Store these instructions for use in a safe place and make them accessible to all users.
- Make sure you include these instructions for use when handing over the device to third parties.

Electric shock

Warning!

- As with every electrical device, use this massage chair carefully and cautiously in order to avoid hazards resulting from electric shock.
- Use the device only with the mains voltage specified on the device.
- Never use the device if it or any of its accessories are visibly damaged.
- Do not use the device if it was accidentally splashed with water.
- Do not use the device during a storm.

- In the event of defects or operational faults, switch off the device immediately and pull the mains plug out of the socket.
- Never pull on the mains cable to remove the mains plug from the socket.
- Do not clamp or pinch the mains cable.
- Keep the mains cable away from hot surfaces.
- Ensure that the device, switch, power supply unit and cable do not come into contact with water, steam or other liquids.
- Use/operate the device only with dry hands.
- Do not subject the device to any impacts or drop the device.

Repairs

Warning!

- Do not open the device.
- Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. For repairs, please contact Customer Service or an authorised retailer. Failure to do so invalidates the warranty.
- If the mains plug or mains cable of the device is damaged, it must be replaced by the manufacturer or at an authorised customer service point.
- The zip fastening of the massage chair must not be opened. It is only used for technical production reasons.
- Do not insert needles or other sharp objects into the chair.

Risk of fire

Warning!

There may be a risk of fire if the device is used improperly or these instructions for use are ignored!

- Ensure that the vents are free at all times. Always keep the vents free from dust, hair etc.
- Do not use the device together with other electrical devices (e.g. electric blanket).
- Never use the device near petrol or other highly flammable substances.

Operation

Warning!

- The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.
- This device is not intended for use by people (including children) with restricted physical, sensory or mental skills or a lack of experience and/or a lack of knowledge, unless they are supervised by a person who has responsibility for their safety or are instructed by this person on how to use the device.
- Do not touch the massage mechanism or the rollers with your hands during operation; there is a danger of crushing.
- Do not insert any objects in the device apertures or in the rotating parts. Ensure that the moving parts can move freely at all times.
- Children should be supervised around the device to ensure they do not play with it.
- Always switch off the device at the main switch after use.
- Ensure that the mains cable does not pose a trip hazard.
- Never stand on the device or place any objects on the device.
- Never massage areas that exhibit swellings, infections or eczema.
- Do not use the device in the event of inexplicable stomach ache.
- Only one person may use the massage chair at any given time.
- Never load the device with more than 120 kg.

Caution!

- Do not expose the device to direct sunlight.
- Do not use the device at temperatures above 40°C.
- Check the function of the device if it has not been used for a longer period of time.


3. Getting to know your device

You can enjoy your own massage at home with the Beurer MC4000 massage chair. The massage chair offers different types of massage to give you a soothing and revitalising massage.

There are three massage programmes that focus on different areas. In addition, the massage chair also features an air pressure massage in the seat and calf area.

4. Unit description

External features

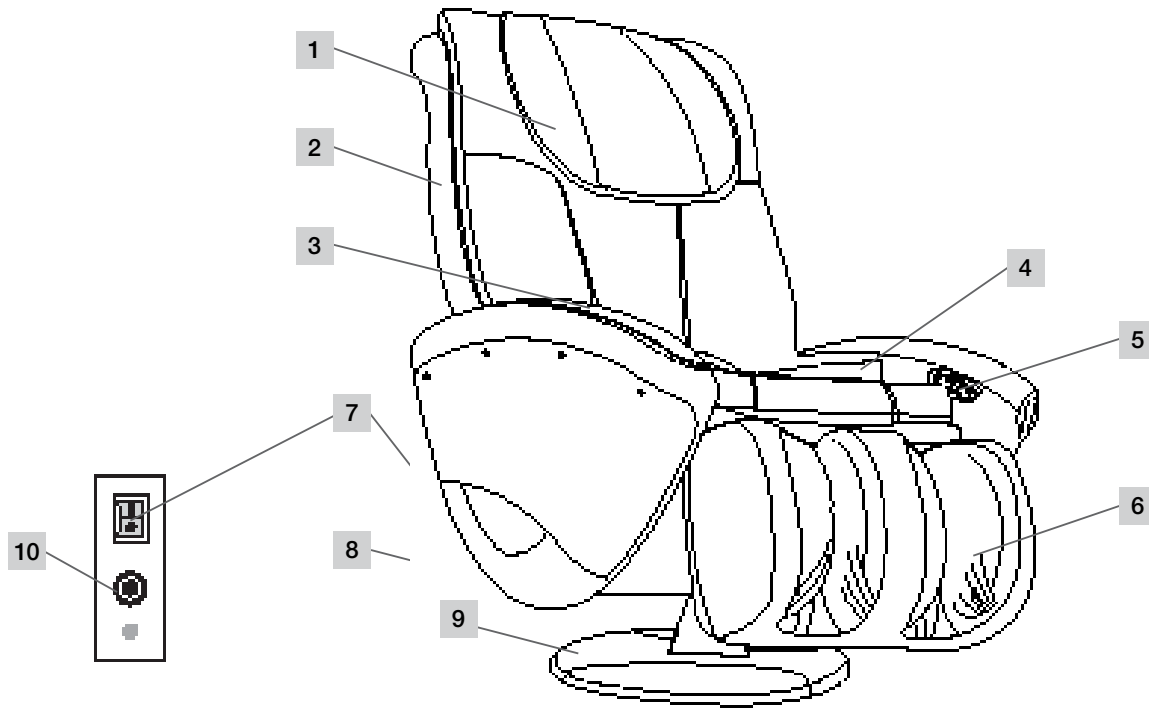1 Neck cushion

2 Back rest

3 Arm rest

4 Seat cushion

5 Controls

6 Leg rest


7 Main switch

8 Locking mechanism

9 Pedestal


10 Cutout switch

Internal features

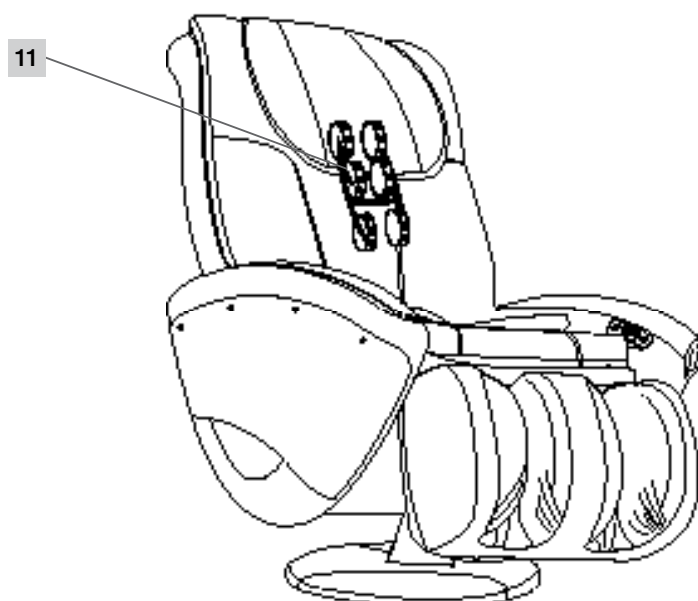11 Massage rollers

Controls

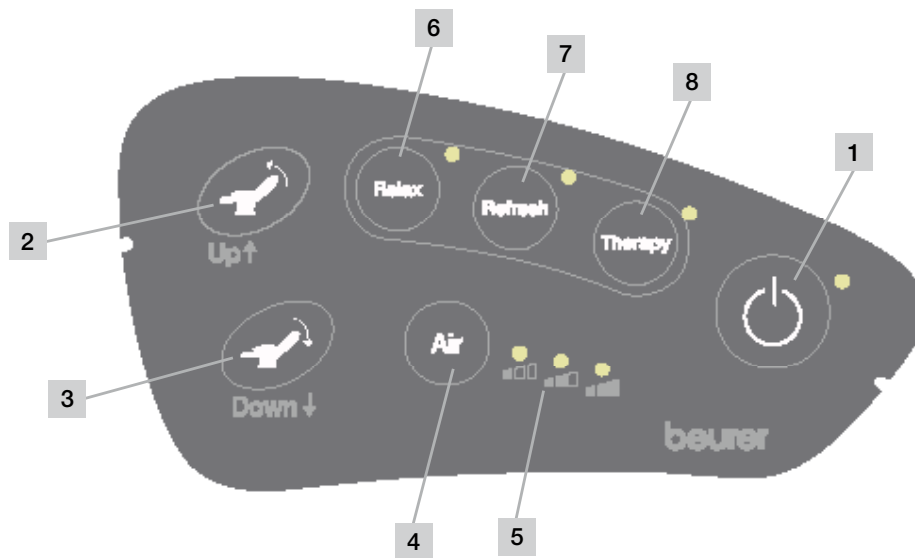- | | |
|-----------------|---|
| 1 On/off button | 5 Intensity of seat and calf air pressure massage |
| 2 Up button | 6 Relax button |
| 3 Down button | 7 Refresh button |
| 4 Air button | 8 Therapy button |

Controlling button functions

On/off button

The controls are activated using the on/off button.
The display is illuminated when the massage chair is in operation.
In standby mode, the display flashes.

Up button

Pressing the up button raises the back rest and lowers the leg rest at the same time.
Acoustic signals sound during the motion.
▷ Upright position/flat position section

Down button

Pressing the down button lowers the back rest and raises the leg rest at the same time.
Acoustic signals sound during the motion.
▷ Upright position/flat position section

Air button

There are four settings for the air pressure massage:
• low • medium • high • off
The currently selected setting is displayed.

Massage programmes

Relax button Refresh button Therapy button

The massage chair features three different automatic massage programmes. These automatic programmes massage your entire back with a combination of kneading, tapping, rolling and Shiatsu massage elements. The air pressure massage can be activated or deactivated separately.

5. Initial use

Warning!

- Keep packaging material away from children! There is a choking hazard.
- Before using the device for the first time, read the safety notes in the Safety section.
 - ▷ Safety, pages 11 to 12

Caution!

- The massage chair could damage your flooring. Place a suitable, non-slip underlay under the massage chair.

Installation location

- Install the massage chair on a sturdy and even surface.
- Before using the device for the first time, remove all obstacles in front of and behind the massage chair so that the back rest can be lowered and the leg rest raised freely.


Unpacking and assembly

The massage chair MC4000 is delivered in a cardboard box.

Caution!


Do not use sharp or pointed tools (Stanley knife) to remove the packaging.

Step 1:

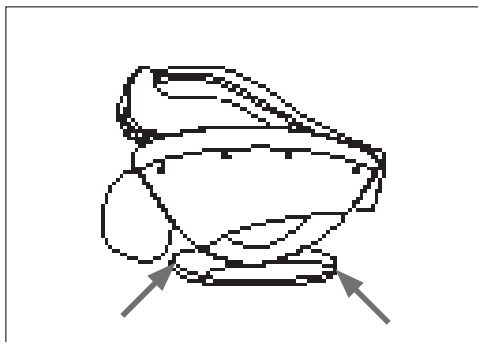
- > Open the cardboard box and hold the massage chair at the designated holding points.
- > Remove the chair from the cardboard box.
Four persons are required for this.
- > Carefully remove all the packaging material.

Step 2:

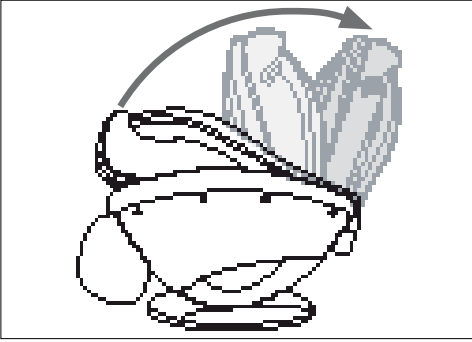
- > Unfold the back rest to the upright position until the locking mechanism audibly engages.
- The back rest now stands at approx. 115°.*

Step 3:

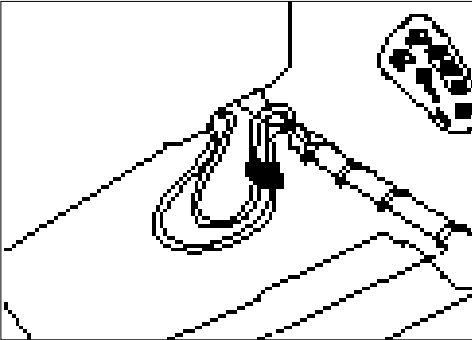
- > Remove the seat cushion.
- > Connect the air hoses under the seat with the air hoses on the seat cushion.
- > Affix using the Velcro strips and ensure that the hoses are not crossed and or kinked.
- > Carefully position the seat cushion.

Step 4:

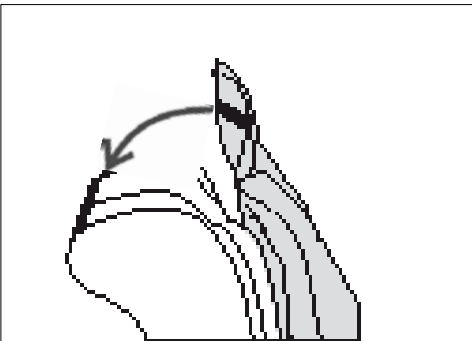- > Affix the back cushion and the neck cushion to the back of the back rest at the top.
- > Affix the Velcro fastening before closing the zip.

Step 5:

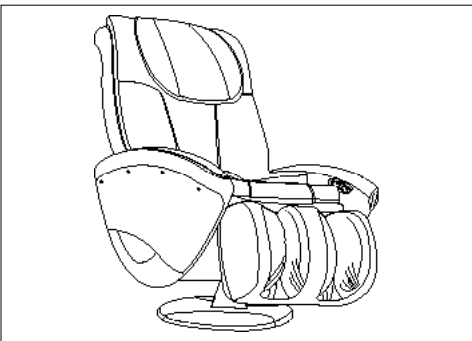- > Plug the mains plug into a socket.
- > Switch on the main switch at the bottom of the rear side.
- > Switch on the massage chair at the  on/off button on the controls.

The massage chair is now ready for operation.

6. Operation

Warning!

- Before using the device for the first time, read the safety notes in the Safety section.
 - ▷ Safety, pages 11 to 12

Caution!

- If you experience irritations of any kind, stop treatment immediately.

Important!

- Never remove the pillow or back cushion during the massage.

Before switching on


- Wear comfortable clothes during the massage.
- Remove all objects from your trouser pockets (e.g. keys, pocket knife, mobile phone).
- Assume a comfortable, relaxed position.
- Always sit in the middle of the seat so that your spine is positioned precisely between the massage rollers.
- The massage should feel pleasant and relaxing at all times. If the massage is painful or uncomfortable, stop the massage or change your sitting position.
- If the device has not been used for a longer period of time, check that the device is functioning correctly before use.

Starting the massage

> Start the massage using the  on/off button. After it is switched on, the massage chair goes into a semi-flat position and begins with the refresh massage programme. The air pressure massage is switched on at low intensity.

- To change to a different massage programme, press the relax or therapy button.
- Press the  air button to alter the intensity of the air pressure for the seat- and calf-area massage or switch it off completely.

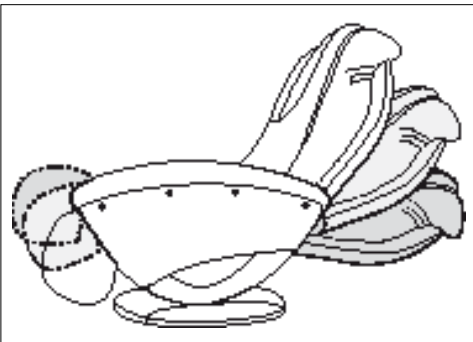Flat position/upright position

The sitting position can be adjusted even when the massage function is switched off.

Putting the massage chair in the flat position

- Press the  down button to put the massage chair into the flat position. This lowers the back rest and raises the leg rest at the same time.
- > Press the  down button once again as soon as the chair reaches the desired position. The motion stops; back rest and leg rest stay in this position.

Putting the massage chair in the upright position


- > Press the  up button to put the massage chair into a higher flat position. This raises the back rest and lowers the leg rest at the same time.

Switching off the massage chair

After the massage programme is finished, the chair automatically enters standby mode.

Caution!


Do not use the device for longer than 20 minutes at a time otherwise the device could overheat. Allow the device to cool for at least 20 minutes before using it again.

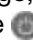During the massage, you can switch off the device at any time using the  on/off button. The massage chair returns to the start position. The operating display flashes.

- > Switch off the massage chair using the main switch beneath the back rest.
- > Pull the mains plug out of the socket.

Additional settings, operating options

Folding mechanism

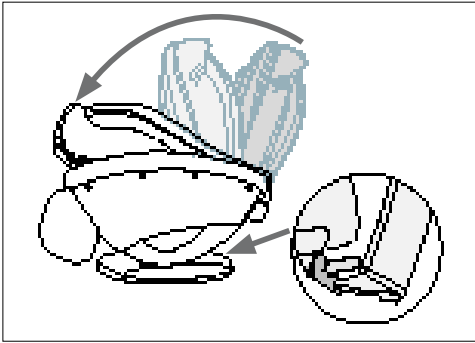
The massage chair can be folded away to save space when stored.

> Switch off the controls at the  on/off button.

The massage chair returns to the start position.

- > Once the chair has reached the start position, switch off the massage chair using the main switch.
- > Pull the mains plug from the socket.
- > Hold the back rest from behind and release the locking mechanism on the underside of the back rest using your foot.
- > Carefully push the back rest forwards on to the seat.
If required, lay the back rest cushion back over the back rest before doing so.

Rotating mechanism

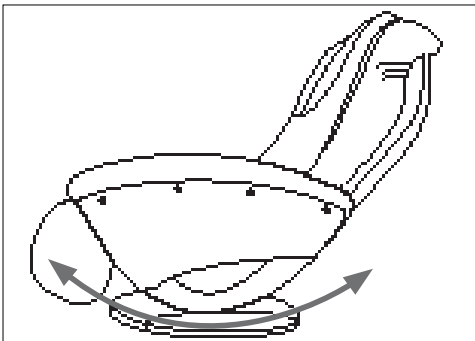The massage chair can be rotated 60° to the left or to the right.

7. Maintenance and storage

Warning!

- The device must be switched off at the main switch after every use and before cleaning.
- Pull the mains plug out of the socket.

- > Wipe the massage chair from time to time with a dry, lint-free or slightly damp cloth.
- > Every two to three months, wipe the chair using a clean cloth with a small amount of care product for artificial leather.
- > Afterwards, dry the chair with a clean, dry cloth.

Artificial leather is sensitive to grease. Clean the massage chair only using care products designated for artificial leather.

Important!

- Clean the device using only the methods specified.
- Under no circumstances may liquid enter the device.
- Do not use any solvent-based cleaning products.
- Do not use aggressive cleaning products.
These can damage the surface.

▷ Folding mechanism, see above

8. What if there are problems?

Problem	Cause	Solution
Grinding noises during operation	These noises are normal and originate from the operating motor or from the massage rollers rubbing against the cushion.	No action is necessary.
The massage cannot be started.	The mains plug is not connected to the socket.	> Plug the mains plug into the socket.
	The massage chair is not switched on at the main switch.	> Switch on the massage chair at the main switch.
	The cutout switch for the massage chair has interrupted the power supply.	> In this case, please contact Customer Service.
Back rest or leg rest cannot be raised or lowered.	Back rest and/or leg rest are overloaded or blocked.	> Remove the objects blocking the device.

9. Technical Data

Product description	Massage chair MC4000
Mains voltage/frequency	220-240 V ~ 50/60 Hz
Energy consumption	150 W
Dimensions (in flat position)	Length 175 cm, width 80 cm, height 67 cm
Weight	72 kg
Maximum load capacity	120 kg

10. Disposal

Observe the local regulations for material disposal. Dispose of the unit in accordance with EC Directive 2002/96/EC – WEEE (Waste Electrical and Electronic Equipment).

If you have any questions, please contact the local authorities responsible for waste disposal.

11. Warranty and service

In the event of a warranty claim or service case, please contact Customer Services in your respective country.

You can find the address in the enclosed address list.

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