# promed 3020

## Promed GmbH

Lindenweg 11

D-82490 Farchant

Tel: +49 (0)8821/9621-0

Fax: +49 (0)8821/9621-21

[info@promed.de](mailto:info@promed.de)

[www.promed.de](http://www.promed.de)

Gebrauchsanleitung | Instruction leaflet | Instruction d'utilisation | Instrucciones de empleo | Istruzioni per l'uso  
Gebruiksaanwijzing | Руководство по эксплуатации | Instrukcja użytkowania | Bruksanvisning | Käyttöohje

# promed 3020

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SE

FI

- A** Mains lead socket
- B** Hand unit lead socket
- C** Quick-action chuck
- D** On/off switch
- E** Rotation speed control/rotation direction control
- F** Digital display
- G** Hand unit holder
- H** Socket for foot pedal lead
- I** On/Off foot pedal, optional

Before replacing the grinding head, the speed control must be set to zero. Failure to observe this instruction may result in significant damage to the handpiece.

1. Turn the quick lock (**C**) to the left.

The chuck will open.

2. Replace the grinding head.

3. Then turn the quick lock (**C**) to the right. The chuck will close and lock the grinding head. To use the device, switch it on again.

## DESCRIPTION OF THE UNIT

Please keep this user manual in a safe place!

Please read the user manual thoroughly before using the appliance for the first time.

Thank you for purchasing the **Promed 3020** manicure and pedicure unit. This unit has been designed and manufactured in accordance with the Medical Device Directive 93/42/EEC. Both qualified specialists and private persons may use the unit after they have read the user manual. This unit is only suitable for use with manicures and gentle pedicures (e.g., removal of thickened skin and calluses, cuticle care, nail care, treatment of thickened nails, etc.). We, as the manufacturer, cannot be made liable in any way for injury or damages to people or objects that arise from failure to comply with this warning. The unit is also suitable for use in the area of nail techniques and nail design as well as the handling of all materials used in this area such as gel and acrylic.

We wish you a lot of enjoyment with your new **Promed 3020** unit. We would like to familiarize you with the unit in the following sections. Please read the user manual thoroughly before using the unit for the first time.

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

## INFORMATION ON THE UNIT / ACCESSORIES

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

Before operating the unit for the first time, please check that all contents have been delivered and that the line voltage is correct. This is done by checking whether the line voltage required for the unit (you can find this information on the specification plate) corresponds to the local power supply (operating voltage of 230 V ~ connected to a correctly installed wall socket with a supply voltage of 230 V ~; operating voltage of 115 V ~ connected to a correctly installed wall socket with a supply voltage of 115 V~).

- The **Promed 3020** unit, in combination with the grinding heads designed for this unit, is ideal for manicures and gentle pedicures, and is also particularly suitable for use with nail techniques (nail design).
- The **Promed 3020** unit uses a high performance motor. The electronic speed control is used to smoothly set the speed of the grinding head. The unit can rotate clockwise or anticlockwise.
- The grinding head holder is able to fixate all commercial grinding heads with standardized dental size shafts of 2.32 to 2.35 mm (0.091 - 0.093 inches). Please use only high-quality grinding heads with a rounded shaft end and a maximal shaft length of 36 mm to guarantee the longest life possible for the grinding head holder! Always make sure that the grinding head does not come in contact with other objects when you lay the handpiece down. There is a risk of the shaft becoming compressed and asymmetrical. Asymmetry in the grinding head causes damage to the handpiece!

## INFORMATION ON THE UNIT / ACCESSORIES

- Grinding head bits are available at your local specialist retailer or **Promed**. ([www.promed.de](http://www.promed.de)).
- To ensure hygienic operation, the grinding heads must be disinfected after use.
- The handpiece may only be disinfected using disinfectant wipes.

Information about suitable disinfectants and products can be obtained from **Promed**.

## OPERATING THE UNIT

Before you switch the unit on, you must set the rotating speed (rotary knob **(E)**) to zero.

The power plug is intended as a separator. Make sure that the power plug is plugged into an easily reachable wall socket.

**A.** Before you start using the device, check if the voltage indicated on the device corresponds to the local mains voltage. Then insert the power plug into the power outlet and the other end of the power cable into the matching socket **(A)** on the left side of the device.

**B.** Insert the plug of the handpiece into the socket **(B)** on the right side of the device. The socket **(H)** next to it is intended for the connection of the optional foot pedal.

**C.** Select a suitable grinding head. Open the quick lock **(C)** by turning it counter-clockwise. Insert the grinding head on the front of the handpiece into the holder for the grinding head. Lock and

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

## OPERATING THE UNIT

secure the grinding head by turning the quick lock clockwise. Never use force to remove or insert a grinding head! Grinding heads with a shaft diameter of 2.32 mm to 2.35 mm fit precisely into the holder for the grinding head without application of excessive force. Please switch off the device when removing a grinding head. Turn the quick lock counter-clockwise and take out the grinding head.

**D.** To switch on the device, press the toggle switch labelled with "Power" into the "ON" position. The device is switched on when the display shows 00000. To switch off, press the same switch into the "OFF" position.

### E. Control of the rotation direction and speed:

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

## OPERATING THE UNIT

### 1. Selection of the rotation direction:

You make this selection by turning the rotary knob **(E)** to the left or right. The rotation direction will also be shown on the display.

Counter-clockwise: the device runs backwards and the display shows 

Clockwise: the device runs forwards and the display shows  .

### 2. Selection of speed:

The speed is regulated by turning the rotary knob **(E)** to the left or to the right. The more you turn, the higher the rotation speed.

A hand symbol will simultaneously appear on the display next to the speed and rotation direction.

### 3. Use of the foot pedal:

Connect the foot pedal **(I)** to the socket **(G)**. A foot symbol will appear on the display. Set the rotation direction and speed by turning again the rotary knob **(E)**. The device is now in the standby mode. Switch on the device by pressing the foot pedal. The grinding head will begin to rotate in the direction and at the speed selected by you. By pressing again the foot pedal, you can switch again into the standby mode. To switch off the device completely, move the power switch **(D)** into the "OFF" position.

#### **Note:**

You may use only **Promed** foot pedals.

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

## OPERATING THE UNIT

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

### **In cases of malfunctioning:**

- Check that the power supply is correct.
- Check whether the handpiece is connected correctly to the unit.
- Check if the device is in the "ON" position (the display should show 00000).
- Check if the foot pedal is correctly connected to the device (the display must show a foot symbol).
- If the problems have still not been solved after following these steps, please consult your service center immediately.

## SPECIAL QUALITY FEATURES

### **ASC: Automatic speed control (digital load compensation control):**

The electronic system registers the pressure exerted on the attachment during filing and compensates for the reduction in rotation speed, so that you always have the maximum power at your disposal.

### **OPC: Overload protection control**

If the grinding head is blocked or stopped for too long, the electronics registers this event and switches the device to a lower, safer operation level. However, you should not keep the device in this state over a long period of time because it could damage the electronics.

### **Maintenance-free, brushless high-performance micromotor:**

The brushes in conventional motors are subject to constant wear and tear until contact finally breaks down altogether, necessitating repairs. Brushless motors are maintenance-free and also offer approx. 20% better performance than motors of a comparable size that use brushes.

### **All-in-one control knob:**

For ease of use, you can regulate both direction and speed of rotation with just the one control knob.

### **LED display:**

The large, easy to read display provides you at a glance with all the information you need on rotation speed, direction of rotation and manual/foot pedal control.

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

## MAINTENANCE TIPS

- A.** Keep the inner parts of the handpiece as free as possible from dust and grinding residue.
- B.** Remove dirt from the outer surface of the handpiece with a cloth moistened with alcohol. Prevent moisture from entering the unit.
- C.** Only use grinding heads with a dental standard size of 2.32 - 2.35 mm (0.091 - 0.093 inches) and avoid using rusted or damaged grinding heads. Their false usage may lead to excessive wear and tear in the ball bearings. Adhere to the safe speed and the safety measures stipulated by the manufacturer.
- D.** Do not expose the handpiece to strong knocks as this may damage the ball bearings. If the handpiece should fall down, check whether the grinding head is bent. Pay closer attention to noise and/or heat generation. Send the handpiece immediately in for repair should it be damaged.

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

## SAFETY INSTRUCTIONS

### **Please pay attention to the following before first time operation:**

Never place the unit in direct contact with fire, gas, oxygen or hot objects such as hot plates.

### **Please observe the following during operation:**

- A.** The handpiece may only be connected to the unit included in the delivery contents.
- B.** Grinding heads to be used:
  - Standardized dental size: 2.32 - 2.35 mm / 0.091 - 0.093 inches.
  - Only use centered grinding heads with a rounded shaft end.
  - The use of rusted or bent grinding heads is not permitted.
- C.** Use the speed controller carefully.
- D.** Switch the unit off before you change the direction of rotation.
- E.** Before you switch the unit on, you must set the rotating speed (rotary knob **(E)**) to zero.
- F.** Appliance power supply: 115 V~ or 230 V~.

- G.** Switch the unit off before you change the grinding head.
- H.** Take every precaution possible to ensure that the unit does not fall down or get damaged in any other way.
- I.** If problems do arise with the handpiece, take it immediately to get it repaired.
- J.** Do not lubricate or wash the handpiece.
- K.** Do not switch the unit on while the chuck is opened.
- L.** Place the test grinding head in the chuck when storing the handpiece.
- M.** Never remove the dust covers from the switches as we can otherwise no longer guarantee that your unit will function properly.
- N.** Only for home use.
- O.** Scrupulous supervision is required when the device is being used in the vicinity of children.

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

## SAFETY INSTRUCTIONS

### **DANGER!**

- Do not place the unit in direct contact with water or other fluids.
-  Never place or use the unit in a wet or moist environment.
- Do not touch the unit with wet hands.
- Only operate the unit with the specified line voltage.
- It is not permitted to use the unit outdoors.

### **WARNING!**

- Make sure that children do not play with the device.
- Never leave the unit unattended when children or inexperienced persons are in its vicinity.
- Only use the unit for the applications described in this manual.
- Only use attachments and accessories from the manufacturer.
- Never use the unit after the electrical lead has been damaged,

when the unit has fallen down or when it has come in contact with water.

- This unit is not designed to be used by people (including children) who have limited physical, sensory, or intellectual abilities or who do not have sufficient experience and/or sufficient knowledge. This does not apply if these people are supervised by a person responsible for their safety or receive instructions from them on how to use this unit.
- Never let children play with the packaging material; there is a risk of suffocation.
- An electrical appliance is not a toy. Children do not recognize the risks that could develop when using electrical appliances. Therefore, operate and store the unit out of the reach of children.
- Do not place or use the unit on hot surfaces or near naked gas flames to ensure that the casing does not melt.

DE

US

FR

IT

ES

NL

RU


PL

SE

FI

## SAFETY INSTRUCTIONS

- Never let the unit come into contact with water. Use a moist cloth to clean the unit. Do not use solvents.
- Never carry or pull the unit by the power cord. Use the plug to pull the power cord out of the wall socket. Do not pull on the power cord itself.
- Do not pull the power cord over sharp edges, wedge it in between objects, hang it down or position or hang it over hot objects (e.g., hot plates) or naked flames. Protect the power cord from heat or oil.
- Only use extension cords when you have satisfied yourself beforehand that they are not defect or damaged.
- If the unit's power cord should become damaged, always replace it with a new one. (Available in electrical specialist stores).
- Do not store the unit near a wash basin or bath tub as there is the danger that the unit may fall or be pulled into the wash basin or the bath tub.
- To avoid all risks of injury, blockage or damage to the unit, make sure that the unit does not come into contact with clothing, cords, shoelaces, bristles, wires, threads, etc. when in operation.

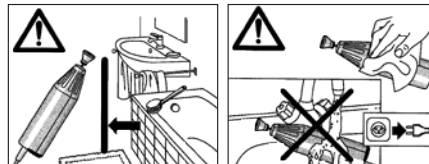DE

US

FR

IT

ES

NL

RU

PL

SE

FI

## SERVICE

Your customers will definitely be expecting you to deliver a professional result. This is only possible with an appliance that is neither defect nor damaged. Therefore, take advantage of our service and have your appliance regularly examined by our specialists. Send in your appliance to Customer Service every 12 months at the latest. We would be happy to provide you with a quote for the costs charged for appliance servicing.

## DISPOSAL

Power tools, accessories, and packaging should be recycled in an environmentally suitable manner.

### Only for EU countries:

Do not dispose of power tools with the domestic waste! In accordance with the European Directive 2002/96/EEC on waste electrical and electronic equipment and its implementation into national law, nonfunctioning power tools must be collected separately and recycled in an environmentally suitable manner.

This symbol indicates that disposal of this product in the domestic waste is not permitted within the EU. Waste devices contain valuable recyclable material that should be recycled. Recycling also prevents uncontrolled waste disposal from damaging the environment and human health. Therefore, please dispose of all waste devices in the appropriate collection systems or send the device to the place where you bought the device for disposal. They will then recycle the device.

DE

US

FR

IT

ES

NL

RU

PL

SE

FI

## WARRANTY

This unit has been manufactured with the utmost of care and has been examined in detail before leaving the factory. We therefore provide a warranty of 24 months after the date of purchase when presented with the warranty card issued with this appliance, in accordance with the following conditions.

- Within the warranty period, we will replace all defect parts of the appliance free of charge, including labor costs for repairs under warranty, in cases of verifiable faults in material or manufacture that have occurred with correct usage and have been identified during the warranty period.
- This warranty is void when the unit has been improperly used (e.g., connection to unsuitable sources of electricity, breakage), the unit has been tampered with (e.g., opening the unit's casing), and when spare parts have been used that have not been approved of by **Promed**. Parts subject to wear are exempt from the warranty. Appliance motors and mobile parts are not included in the warranty.
- The warranty period begins on the date of purchase. Warranty claims do not affect the warranty duration. Warranty claims must be lodged within the warranty period. Claims arising after cessation of the warranty period cannot be considered.
- The warranty only comes into effect within the scope of these warranty conditions when the date of purchase has been confirmed by the stamp/ signature of the dealer on the warranty card.
- In cases of warranty or repair claims, please post the complete appliance with the completed warranty card to the customer service point responsible for you. The warranty card can be found at the end of this user manual.
- We reserve the right to make technical and optical changes as well as changes to the equipment!

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### Symbol for protection class II appliances

The symbol corresponds to symbol 10, table D1, annex D of the EN 60601-1:1990 Standard. The appliance is double insulated as a protective measure. This prevents conductive parts of the appliance from becoming impaired, something that could lead to malfunctioning of the basic insulation voltage.

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### Symbol for type B

The symbol corresponds to symbol 1, table D2, annex D of the EN 60601-1:1990 Standard. Appliances that correspond to the requirements of this norm offer a higher level of protection against electric shocks, especially when taking into account the permissible patient leakage currents and the safety of the user.

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### Symbol for WARNING! Indicates sections of the user manual that are relevant for safety!

The symbol corresponds to symbol 14, table D1, annex D of the EN 60601-1:1990 Standard. Instructions for an appliance, system or accessory that contain important information for the user that are relevant for safety.

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### Pay attention to the symbol for the user manual!

The symbol corresponds to the D-M018 symbol in the DIN 4844-2 Standard. It is imperative that the manual be read before using the unit.

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### US | CE-Directives

The appliance is equipped with radio interference suppression in accordance with the EC Directive 87/308/EEC. It is designed for operation with a line voltage of 230/115 V~, 50/60 Hz. CE according to the EC Directive 2006/95/EEC in accordance with the EC Directive 2004/108/EEC and EC Directives 2002/95/EEC, 1907/2006/EEC and 2005/69/EEC.

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The requirements of the EC Directive 93/42/EEC for class IIa medical products has been fulfilled according to the Standard EN 60601-1:1988+A1:1991+A2:1995; EN 60601-1:2003. The CE marking for the appliance refers to the EC Directive 93/42/EEC. Appliance classification: Class II.  
In accordance with EMC testing in  
EN 60601-1-2:1994; EN 55011 KLB:1998;  
EN 61000-4-2:1996; EN 61000-4-3:1997;  
electromagnetic disturbances are neither emitted from the appliance nor affect the appliance.

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